

# ENDO



## ENDO 550S

Surgical Table Hydraulic

<https://endo.id>

# ENDO 550S

## Surgical Table Hydraulic

ENDO Surgical Table Hydraulic ENDO 550S are designed to accommodate various type of patient positioning required for demanding surgical procedures with ultimate care to the patient's and user's safety. Surgical Tables are constructed with stainless steel and high quality materials to provide durability and long operating life.

The tabletop consists of four HPL-Compact Laminant sections: head, backrest, seat and leg; each with removable, easy to clean, antibacterial, and waterproof polyurethane pads for patient comfort.

The radiographic tabletop\* made of HPL - Compact Laminant can be easily installed and removed for X-ray procedures.

The durable cast iron table base is covered with a stainless steel. The "C" shape design of the ENDO 550S allows easy access to the tabletop in any position. Special brakes and antistatic casters allow ease of mobility. Stainless steel side rails accept a variety of available accessories.

ENDO 550S Surgical Tables are reliable products having both quality and competitive pricing it is an excellent choice.



### ENDO 550S Surgical Table, Hydraulic

The tabletop can be raised and lowered to the desired height with a hydraulic lift in a quiet and nonvibrating manner by using a foot pedal. Downward speed is adjustable. Other specific table settings are performed manually. The tabletop can be rotated for 355°.

(\* ) Optional

## Technical Details


- Height adjustment:  
770 - 1000 mm
- Height adjustment (w/mattress):  
800 - 1030 mm
- Tabletop (w/o accessories):  
450 x 2020 mm
- Table base:  
550 x 800 mm Head section: 450 x 310 mm  
Backrest section: 450 x 550 mm
- Seat section: 450 x 445 mm
- Right leg section: 215 x 680 mm
- Left leg section: 215 x 680 mm

## Packing Details

- Net weight: 220 kg
- Gross weight: 265 kg
- Packing dimensions: 1000 x 1680 x 1100 mm

## Positions

- Head section: Upward 15°, downward 90°,  
Extendable 5 cm, Removable
- Lateral: Right 20°, Left 20°
- Trendelenburg: 25°
- Reverse Trendelenburg: 25°
- Backrest section:  
Downward 15°,  
Upward 75° (w/mattress),  
Upward 90° (w/o mattress),  
Upward 20° (w/radiographic tabletop and  
mattress)
- Reflex: 105°
- Flex: 195°
- Leg section: Downward 90°, Upward 15°  
Spreading 90°, Removable, Independent  
movement
- Renal surgery: 165°
- Tabletop: Rotates 355° (550S)

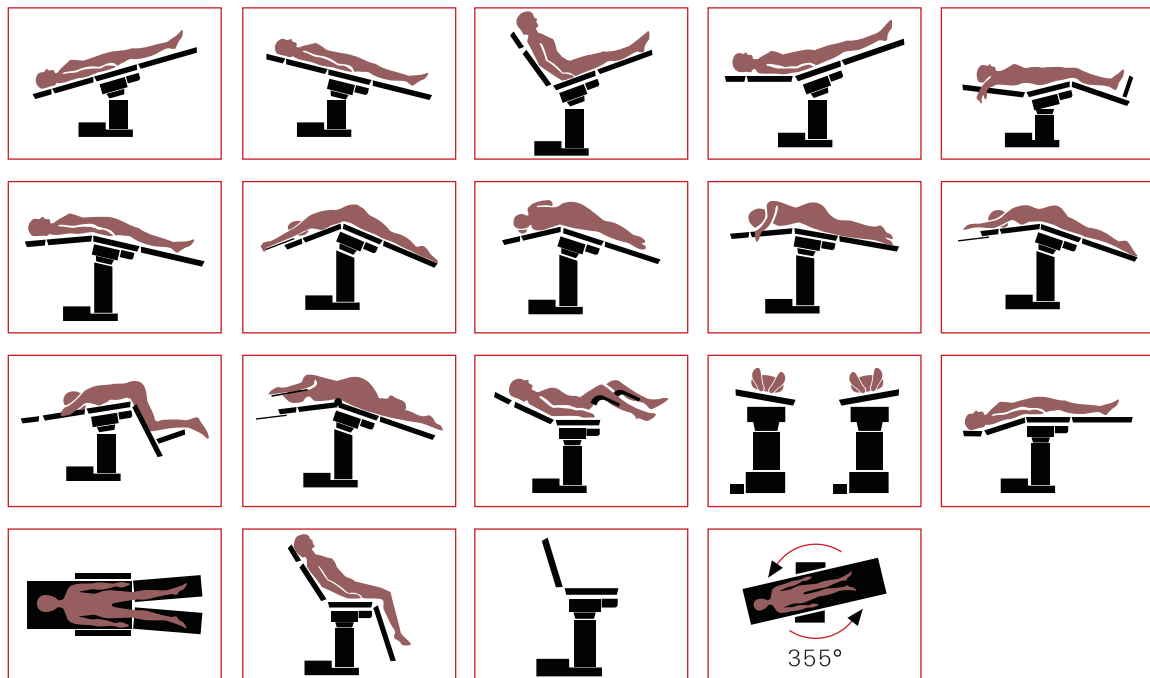

# Accessories For ENDO 2000E, ENDO 550S Surgical Tables

## Standard accessories

- Anaesthetist's frame (1 ea)
- I.V. stand (1 ea)
- Shoulder supports (2 ea)  
(Optional accessories for 2000E)
- Body supports (2 ea)  
(Optional accessories for 2000E)
- Surgical table arm boards (2 ea)
- Knee crutches (2 ea)  
(Optional accessory for 2000E, 550S)
- Compact stand for surgical table accessories (1 ea)  
(Optional accessory for 2000E, 550S)
- Body strap (1 ea)
- Mattress (1 set)



• Anaesthetist's frame  
Ref: M3S 000 0 0



• I.V. stand  
Ref: M3T 000 0 0



• Shoulder supports  
Ref: M3P 000 0 0 (Ref of 1 ea)  
(Optional accessories for ENDO 2000E)



• Body supports  
Ref: M3R 000 0 0 (Ref of 1 ea)  
(Optional accessories for ENDO 2000E)



• Surgical table arm boards  
Ref: M3M 000 0 0 (Ref of 1 ea)



• Knee crutches  
Ref: M9N 000 0 0 (Ref of 1 ea)  
(Optional accessory for ENDO 2000E,  
ENDO 550S)



• Compact stand for surgical  
table accessories  
Ref: M3E 100 5 0  
(Optional accessory for ENDO 2000E,  
ENDO 550S)



• Body strap (1 ea)  
Ref: M9E 000 0 0



• Mattress  
Ref: M2C 800 5 0 (2000E)



• Polyurethane mattress  
Ref: M2C 903 5 0  
Ref: M2D 901 5 0 (550)

## Accessories

- Foot support
- Knee crutches
- Arthroscopic knee crutch
- Knee rest
- Proctology attachment
- Radiographic tabletop assembly kit
- Radiographic tabletop kit (550)
- X-ray cassette holder
- Wrist grip
- Arm strap
- Arm/hand surgical section
- Waste tray
- Leg rest - Lithotomy position
- Head section, high adjustable
- Arm support
- Horseshoe headrest
- Skull clamp
- Crossbar adaptor for skull clamp
- Stand for skull clamp and horseshoe headrest
- Arthroscopic shoulder support
- Lithotomy leg positioning boat
- Spinal surgery frame
- Carbon fiber leg section (2000E)
- Compact laminate leg section (2000E)



• **Foot support**  
Ref: M3J 000 0 0 (Ref of 1 ea)



• **Arthroscopic knee crutch**  
Ref: M2D 309 5 0



• **Knee rest**  
Ref: M9A 001 5 0 (200E)  
Ref: M30 000 0 0 (550)



• **Proctology attachment**  
Ref: M3I 000 0 0



• **Radiographic tabletop assembly kit**  
Ref: M2C 400 5 0



• **X-ray cassette holder**  
Ref: M9D 000 0 0



• **Wrist grip**  
Ref: M2D 333 5 0



• **Proctology attachment**  
Ref: M3I 000 0 0



• **Waste tray**  
Ref: M9H 010 0 0 (2000E)  
Ref: M3Z 000 0 (550)



• **Leg rest - Lithotomy position**  
Ref: M2D 600 5 0 (Ref of 1 ea)



• **Head section, high adjustable**  
Ref: M9B 200 5 0  
Ref: M9B 000 0 0 (550)



• **Head section, high adjustable**  
Ref: M9B 250 0 0 (2000E)



• Arm support  
Ref: M9M 000 0 0 (Ref of 1 ea)



• Horseshoe headrest  
Ref: M3F 000 0 0



• Skull clamp  
Ref: M9A 000 0 0



• Crossbar adaptor for skull clamp  
Ref: M9A 001 5 0



• Stand for skull clamp and horseshoe headrest  
Ref : M9K 000 0 0



• Arthroscopic shoulder support  
Ref: M9S 150 0 0



• Lithotomy leg positioning boat  
Ref: AU-OTJ01-001 (Ref of 1 pair)



• Spinal surgery frame  
Ref: M9P 200 0 0




• Carbon fiber leg section  
Ref: M9S 050 0 0 (2000E)



• Compact laminate leg section  
Ref: M9S 060 0 0 (2000E)

# ENDO

 Jalan Raya Menganti 14 Kedurus, Surabaya 60223  
Jawa Timur, Indonesia

 62-31-7673636

 62-31-7673737

 0800 . 177 . ENDO (3636)

 info@endo.id

 <https://endo.id>