

ENDO



ENDO 2000E Surgical Table Electro Hydraulic

<https://endo.id>

ENDO 2000E

Surgical Table Electro Hydraulic

ENDO Surgical Table Electro Hydraulic ENDO 2000 is designed to accommodate different patient positions required for demanding surgical procedures, with an ultimate focus on patient and user safety. ENDO 2000E Surgical Table is constructed with stainless steel and high quality materials to provide durability and long operating life.

Surgical positions of ENDO 2000E are adjusted by an electrohydraulic system. ENDO 2000E Surgical Table has up/down, sliding, Trendelenburg/reverse, right/left lateral, renal surgery, Trendelenburg and back section upward combination, and back section tabletop positions. Positions are adjusted with remote control, in a quiet and nonvibrating manner at constant speed.

ENDO 2000E Surgical Table is equipped with battery recharge system. Movements of the table are fed from battery circuit. Battery power of the table has a 30 minutes operating time in cases of mains failure providing approximately 30 complete position changes.

The tabletop consists of four sections: head, backrest, seat and leg, each with removable, easy to clean, anti-bacterial, and waterproof 80 mm mattress for patient comfort.

Easy to install radiographic tabletop kit* enables using ENDO 2000E for X-ray procedures. Sliding motion of the tabletop enables ENDO 2000E Surgical Table compatible for use with C-Arm radiological equipment.



The durable cast iron table base and electrohydraulic piston is covered with stainless steel. The "T" shape design of the ENDO 2000E allows easy access to the tabletop in any position. Special brakes and anti-static casters allow ease of mobility. Stainless steel side rails accept a variety of available accessories.

ENDO 2000E Surgical Table is a state-of-the-art product and with both quality and competitive pricing it is an excellent choice for today's demanding surgery theaters.

(*) Optional

Technical Details


- Height adjustment: 725 - 1055 mm
- Height adjustment (w/mattress): 805 - 1135 mm
- Tabletop: 500 x 1880 mm
- Table base: 370 x 830 mm
- Head section: 500 x 240 mm
- Backrest section: 500 x 490 mm
- Seat section: 500 x 510 mm
- Right leg section: 235 x 555 mm
- Left leg section: 235 x 555 mm
- Voltage: 230V, 50Hz
- Remote control voltage: DC 5V
- Motor inputs: DC 24V, 10A Battery: 2 each, DC 12V, 7Ah, maintenance free, lead acid battery
- Battery recharge voltage: DC 27.2V
- Battery recharge time: 18h (full charge)
- Low Battery Level Alarm: Sound alarm

Positions

- Head section: Upward 15°, Downward 90°, Extendable 5 cm, Removable
- Backrest section: Upward 75°, Downward 45°, Upward 45° (w/radiographic tabletop and mattress)
- Leg section: Upward 15°, Downward 90°, Spreading 90°, Removable, Independent movement
- Lateral: Right 20°, Left 20°Trendelenburg: 25°
- Reverse Trendelenburg: 25°Lateral nephrectomy: 135°Reflex: 105°
- Flex: 225°
- Tabletop: Rotates 355°Horizontal tabletop sliding:250 mm

Packing Details

- Net weight: 250 kg
- Gross weight: 295 kg
- Packing dimensions: 1000 x 1680 x 1100 mm

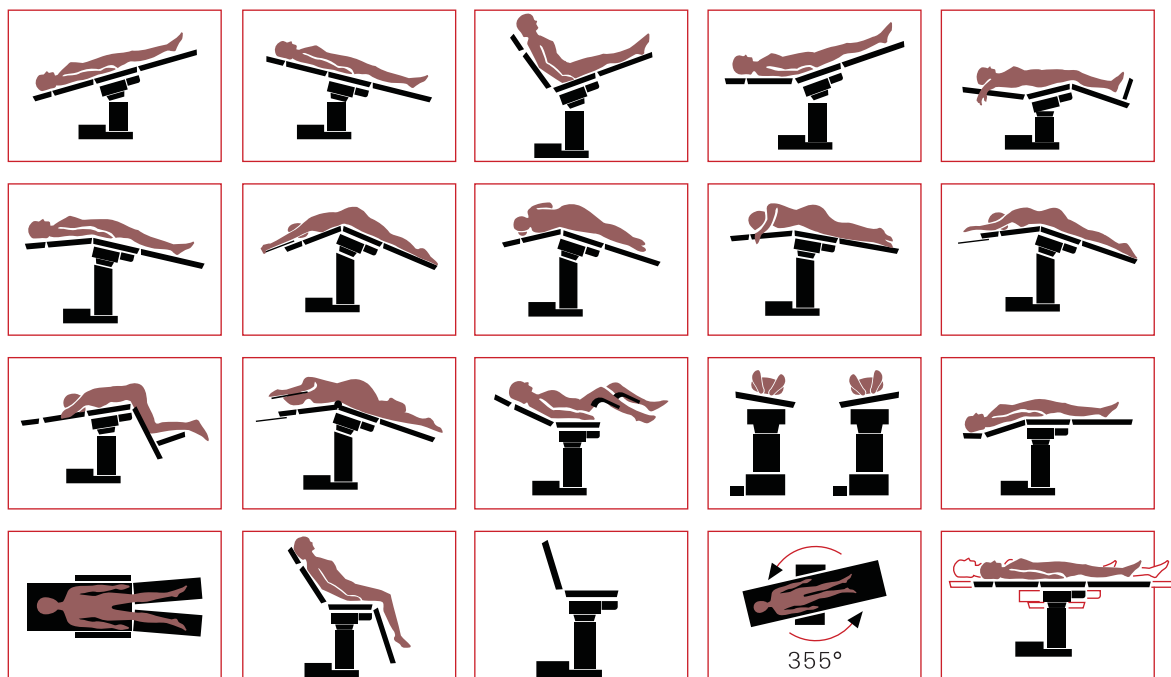

Accessories For ENDO 2000E, ENDO 550S Surgical Tables

Standard accessories

- Anaesthetist's frame (1 ea)
- I.V. stand (1 ea)
- Shoulder supports (2 ea)
(Optional accessories for 2000E)
- Body supports (2 ea)
(Optional accessories for 2000E)
- Surgical table arm boards (2 ea)
- Knee crutches (2 ea)
(Optional accessory for 2000E, 550S)
- Compact stand for surgical table accessories (1 ea)
(Optional accessory for 2000E, 550S)
- Body strap (1 ea)
- Mattress (1 set)



• Anaesthetist's frame
Ref: M3S 000 0 0



• I.V. stand
Ref: M3T 000 0 0



• Shoulder supports
Ref: M3P 000 0 0 (Ref of 1 ea)
(Optional accessories for ENDO 2000E)



• Body supports
Ref: M3R 000 0 0 (Ref of 1 ea)
(Optional accessories for ENDO 2000E)



• Surgical table arm boards
Ref: M3M 000 0 0 (Ref of 1 ea)



• Knee crutches
Ref: M9N 000 0 0 (Ref of 1 ea)
(Optional accessory for ENDO 2000E,
ENDO 550S)



• Compact stand for surgical
table accessories
Ref: M3E 100 5 0
(Optional accessory for ENDO 2000E,
ENDO 550S)



• Body strap (1 ea)
Ref: M9E 000 0 0



• Mattress
Ref: M2C 800 5 0 (2000E)



• Polyurethane mattress
Ref: M2C 903 5 0
Ref: M2D 901 5 0 (550)

Accessories

- Foot support
- Knee crutches
- Arthroscopic knee crutch
- Knee rest
- Proctology attachment
- Radiographic tabletop assembly kit
- Radiographic tabletop kit (550)
- X-ray cassette holder
- Wrist grip
- Arm strap
- Arm/hand surgical section
- Waste tray
- Leg rest - Lithotomy position
- Head section, high adjustable
- Arm support
- Horseshoe headrest
- Skull clamp
- Crossbar adaptor for skull clamp
- Stand for skull clamp and horseshoe headrest
- Arthroscopic shoulder support
- Lithotomy leg positioning boat
- Spinal surgery frame
- Carbon fiber leg section (2000E)
- Compact laminate leg section (2000E)



• **Foot support**
Ref: M3J 000 0 0 (Ref of 1 ea)



• **Arthroscopic knee crutch**
Ref: M2D 309 5 0



• **Knee rest**
Ref: M9A 001 5 0 (200E)
Ref: M30 000 0 0 (550)



• **Proctology attachment**
Ref: M3I 000 0 0



• **Radiographic tabletop assembly kit**
Ref: M2C 400 5 0



• **X-ray cassette holder**
Ref: M9D 000 0 0



• **Wrist grip**
Ref: M2D 333 5 0



• **Proctology attachment**
Ref: M3I 000 0 0



• **Waste tray**
Ref: M9H 010 0 0 (2000E)
Ref: M3Z 000 0 (550)



• **Leg rest - Lithotomy position**
Ref: M2D 600 5 0 (Ref of 1 ea)



• **Head section, high adjustable**
Ref: M9B 200 5 0
Ref: M9B 000 0 0 (550)



• **Head section, high adjustable**
Ref: M9B 250 0 0 (2000E)



• Arm support
Ref: M9M 000 0 0 (Ref of 1 ea)



• Horseshoe headrest
Ref: M3F 000 0 0



• Skull clamp
Ref: M9A 000 0 0



• Crossbar adaptor for skull clamp
Ref: M9A 001 5 0



• Stand for skull clamp and horseshoe headrest
Ref : M9K 000 0 0



• Arthroscopic shoulder support
Ref: M9S 150 0 0



• Lithotomy leg positioning boat
Ref: AU-OTJ01-001 (Ref of 1 pair)



• Spinal surgery frame
Ref: M9P 200 0 0




• Carbon fiber leg section
Ref: M9S 050 0 0 (2000E)



• Compact laminate leg section
Ref: M9S 060 0 0 (2000E)

ENDO

 Jalan Raya Menganti 14 Kedurus, Surabaya 60223
Jawa Timur, Indonesia

 62-31-7673636

 62-31-7673737

 0800 . 177 . ENDO (3636)

 info@endo.id

 <https://endo.id>