


ENDO

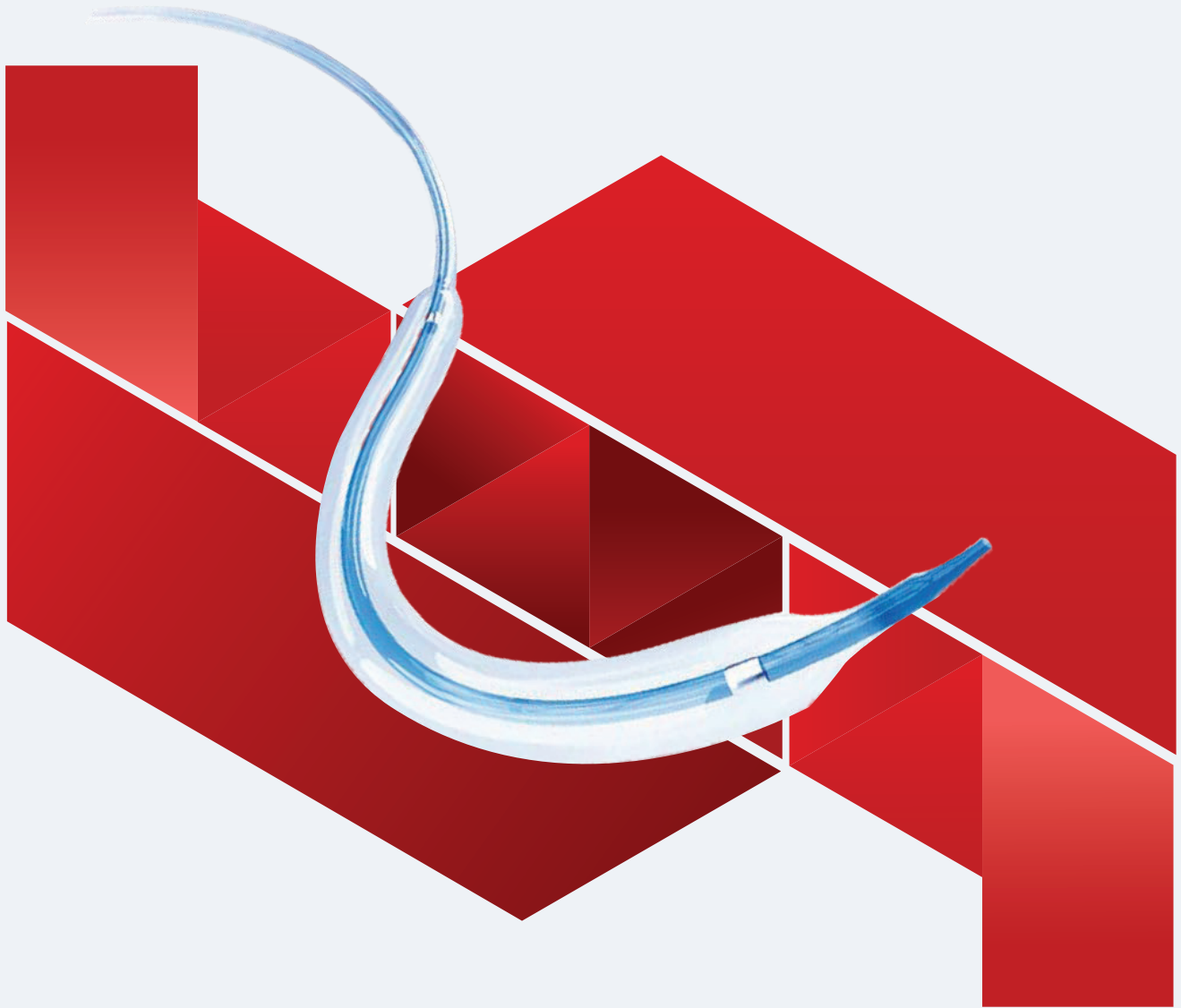

ENDO Coronary Dilatation Catheter (2)

<https://endo.id>

ENDO Coronary Dilatation Catheter (2)

Features:

- Performance for kissing balloon under 6F guiding catheter.
- Patented hydrophilic coating on distal shaft and Hypotube proximal shaft provide good deliverability.
- Three-directional memory folding of the balloon allows of size minimization with great withdrawal capability, and ensures the good performance of repeated dilatation.
- Soft tapered tip provides smooth transition from guidewire to balloon catheter to avoid "fish mouth", and minimize the vessel damage.
- Short shoulder design at both ends of the balloon ensures the precise dilatation, and reduces the damage to normal vessel next to the lesion.

Technical Parameters


Catheter Design	Rapid Exchange
Balloon Compliance	Semi-compliant
Marker	2
Distal Shaft O.D.	2.7F
Proximal Shaft O.D.	1.9F
Tip Profile	0.017"
Minimum Guiding Catheter Compatibility	No less than 5F (0.056")
Guide Wire Compatibility	0.014"
Effective Length	142cm
Nominal Pressure	8 atm (\varnothing 2.25mm-3.0mm) / 6 atm (\varnothing 3.25mm-4.0mm)
Rated Burst Pressure	16 atm (\varnothing 2.25mm-3.75mm) / 14 atm (\varnothing 4.0mm)

Ordering Information


Balloon Diameter (mm)	Balloon Length (mm)						
	8	10	12	15	20	25	30
2.5	ENDO LPCRX 25008	ENDO LPCRX 25010	ENDO LPCRX 25012	ENDO LPCRX 25015	ENDO LPCRX 25020	ENDO LPCRX 25025	ENDO LPCRX 25030
2.75	ENDO LPCRX 27008	ENDO LPCRX 27010	ENDO LPCRX 27012	ENDO LPCRX 27015	ENDO LPCRX 27020	ENDO LPCRX 27025	ENDO LPCRX 27030
3.0	ENDO LPCRX 30008	ENDO LPCRX 30010	ENDO LPCRX 30012	ENDO LPCRX 30015	ENDO LPCRX 30020	ENDO LPCRX 30025	ENDO LPCRX 30030
3.5	ENDO LPCRX 35008	ENDO LPCRX 35010	ENDO LPCRX 35012	ENDO LPCRX 35015	ENDO LPCRX 35020	ENDO LPCRX 35025	ENDO LPCRX 35030
4.0	ENDO LPCRX 40008	ENDO LPCRX 40010	ENDO LPCRX 40012	ENDO LPCRX 40015	ENDO LPCRX 40020	ENDO LPCRX 40025	ENDO LPCRX 40030

ENDO

 Jalan Raya Menganti 14 Kedurus, Surabaya 60223
Jawa Timur, Indonesia

 +62-31-767-3636

 info@endo.id

 +62-31-767-3737

 <https://endo.id>

 0800 . 177 . ENDO (3636)