

ENDO



ENDO DRX

Digital Radiographic X-ray Equipment

<https://endo.id>

ENDO DRX

Digital Radiographic X-ray Equipment

Large size touch user interface Easy exposure control

Graphical anatomy program

The system offers an easy-to-understand visual logic for intuitive exposure parameter setting.

3D positioning guidance

Creating a 3D positioning image to help the patient for fast and accurate positioning guidance.

Self-adaptive screen

Screen auto rotation between table and bucky stand, offering easy and convenient operation.

Data synchronization

Exposure parameters synchronization between touch screen interface and workstation, speeding up your workflow in daily routine.

Emergency mode

Allowing instant exposure for fast examination, saving valuable time for emergency workflow.

DAP monitoring Option

Real time DAP monitoring is available, measuring actual dose level and displaying on the console.



Intelligent system control makes your operation easier



In-tray charging

Built-in charging port in both bucky stand and bucky table can always speed up workflow with fast charging.



Auto stitching Option

Captured image can be automatically stitched together by the system, giving a long image that extend across the larger areas of the body.




Auto tracking Option

Synchronized vertical travel of X ray tube and bucky unit, providing fast positioning for examination.



Auto exposure control Option

Automatically controlling the amount of total radiation to ensure a better image quality.

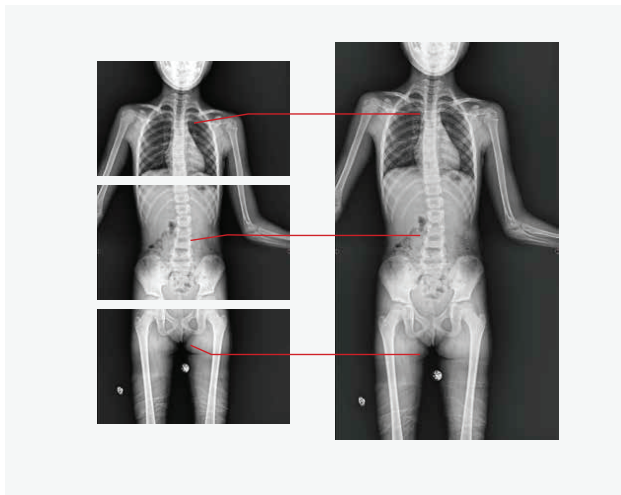

Smooth positioning at complex angle for radiography

Floor-mounted tube support and fixed table

Radiography in both standing position and various orthopedic position.

Flexible X ray tube movement

Free rotation of X ray tube enables different angle of radiography for complex orthopedic position, allowing easy positioning for patient on wheelchair and stretcher.

4-way, floating bucky table


Long table floating range makes easy radiography for every desired anatomy of human body.

Long vertical travel of X ray tube

The X ray tube travels in long vertical distance. 35cm minimum distance between X ray tube and floor make it possible for ankle radiography in weight bearing position.

Flexibility of grids

Equipped with two removable grids for bucky stand (focal distance=1.8m) and bucky table (focal distance=1.0m), dedicating for chest radiography and table radiography.

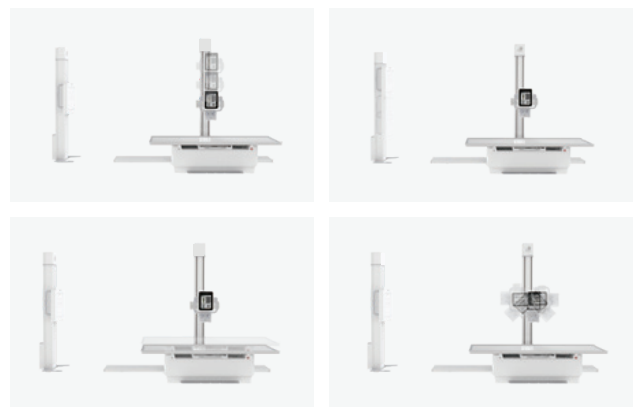

Rich options for clinical application

Multiple bucky table options

Three available options of bucky table: mobile type, fixed type and lifting type. The lifting table is able to lower to a minimum of 50cm, maximizing your focus on patient care.



High inverter frequency generator

- Self-developed high voltage generator provides high quality and stable X ray.
- Unique mobile design maximize your spacial layout.



High performance X ray tube

- X ray tube with high heat capacity and high speed anode rotation.
- Extending your workflow with stable image quality in a variety of clinical application.



High end panel detector

- High spacial resolution imaging detector brings crystal clear image with delicate details.
- Flexible configuration and option: mobile or wireless, single or double, it's all up to you.



Ambitious workstation

- Image workstation expand your diagnosis with versatile features: e.g. cardiothoracic ratio measurement, position measurement.

Selected components contributes to excellent image quality



High efficient image system delivers effortless interaction

Registered exposure programs

The system can register up to hundreds of exposure parameters regard to different anatomies and positions, making it always easy for radiography.

Fast imaging

The system can automatically process the acquired image, bring clear image with speed and improving examination efficiency.

Image post-process

Performing post-processing immediately after acquiring each image based on preset image database, ensuring the best image quality.

Remote data sharing

Seamless interconnection between PACS and image workstation, remote data sharing.


One push on/off switch

Fast boot up at the press of one button, making it extremely simple and convenient for examination.

Barrier-free communication


Radiologist can give instructions to the the patient or adjust the parameters outside of the X ray room through the workstation table.



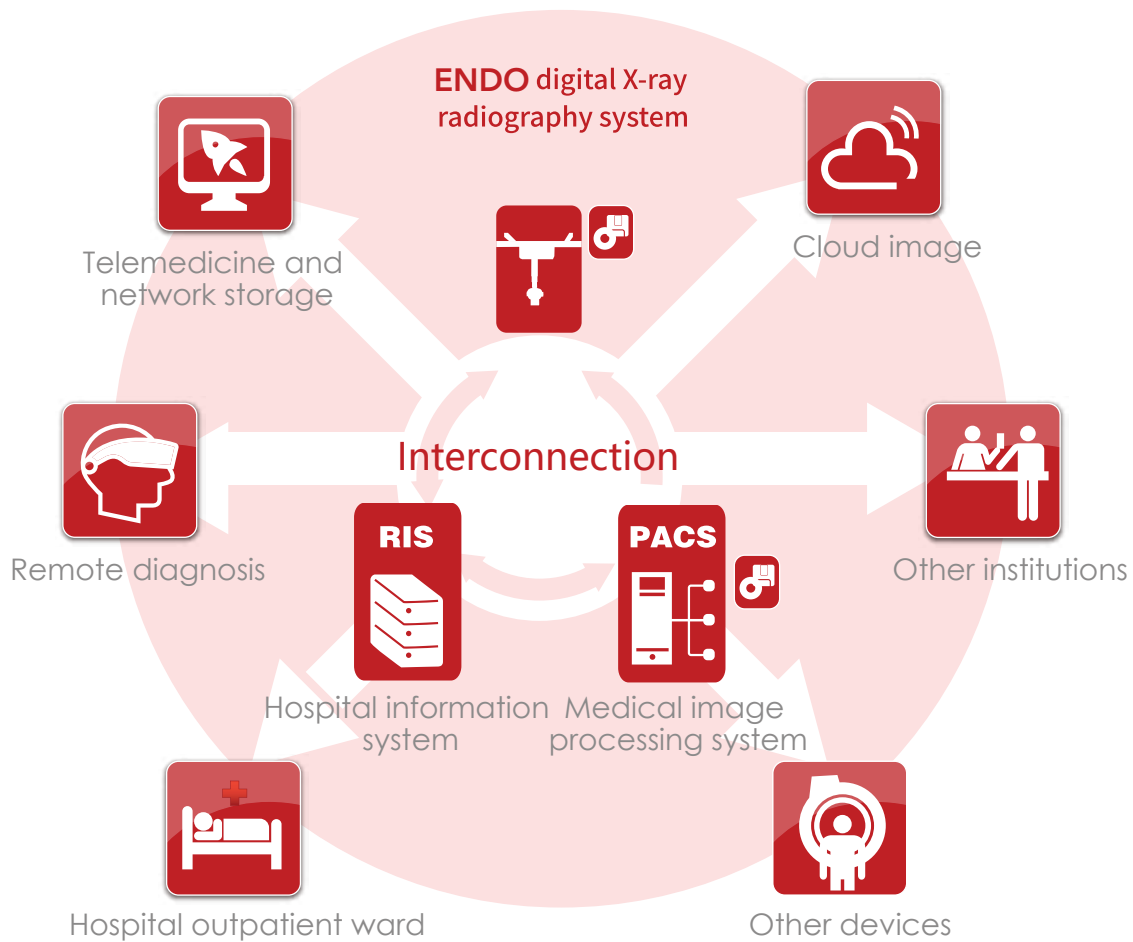


- The system can be maintained remotely, it can perform fault self-diagnosis and display fault codes, quickly and accurately determine the fault status and has self-protection functions.
- It can share information with PACS system, bring benefits for remote diagnosis.


ENDO

 Jalan Raya Menganti 14 Kedurus, Surabaya 60223
Jawa Timur, Indonesia

 62-31-7673636

 62-31-7673737

 0800 . 177 . ENDO (3636)

 info@endo.id

 <https://endo.id>