

ENDO ECG-12

Twelve-channel ECG



ENDO ECG-12


Twelve-channel ECG



Touch Screen



Portable Design

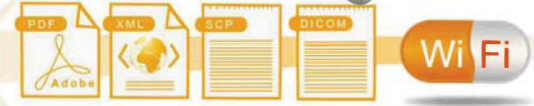

Seamless Connectivity



Hospitals
EMR/HIS/PACs

DICOM/HL7



Clinics
GDT-based

GDT



Smart ECG Viewer



User-friendly Design

- 7" high resolution color touch screen
- Alphanumeric keyboard

Streamlined Workflow

- One-touch operation
- Support barcode scanner to simplify workflow

Accurate Clinical Performances


- Weak signal detection with the advanced 24-bit A/D converter technology
- Anti-noise technology with comprehensive filters and high CMRR over 121dB
- 120s real-time ECG waveform freezing and review to assist diagnosis
- Automatic measurements and interpretation tested with CSE/AHA/MIT databases
- Extend printing once arrhythmia is detected for further rhythm analysis

Flexible Printing Solutions

- Built-in high resolution thermal recorder
- Support external USB printer

Data Management

- Internal storage up to 500 ECGs and can be enlarged by USB flash disk and SD card
- Smart ECG Viewer software for data management on PC (optional)
- Data transmission to PC via series port/LAN/internal Wi-Fi
- Smart ECG Viewer software for data management on PC (optional)
- Report formats: PDF and optional SCP/FDA-XML/DICOM data export

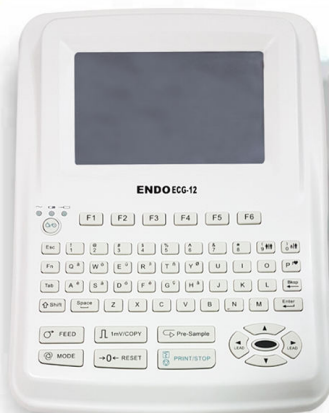

High Resolution Touch Screen with Grid Display



Gender and Age-group Shortcuts



Extensive Communication Interfaces



One-touch Operation

Network your ECG data with the help of ENDO's Smart ECG Viewer in the same LAN. It is convenient to share, review, edit and diagnose ECG reports between different doctors. Meanwhile, with the help of Smart ECG Viewer, ENDO ECG-12 is able to access EMR/PACs/GDT based Systems via the standard protocols: DICOM Worklist/HL7/GDT.

Network ECG data in the LAN


One Smart ECG Viewer connects maximum 10 units ECG machines to manage or edit ECG data.

Hospital System Access

Patient information (worklist) will be transmitted from hospital systems to ECG machines. Meanwhile, reports or final diagnostic results will be uploaded to hospital systems for archiving.

Standard Communication Protocols

HL7/DICOM worklist are supported on Smart ECG Viewer.

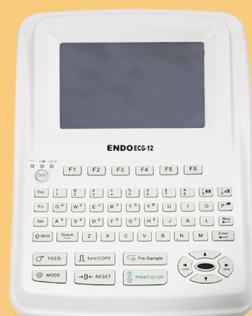

ENDO ECG-12

Twelve-channel ECG

About ENDO

ENDO INDONESIA was founded in 2006 as a company engaged in importation and distribution of Medical Equipment. Today our products include Electromedical Equipment, Reagents, Reusable and Disposable Medical Devices, among many others and accessories. And now we had grown into the manufacturing field of Medical Recording Chart Paper.

ENDO is a company that fulfills legal requirements and receives approval from company's legal entity of Indonesia, with number C-05336 HT.01.01.TH.2006. Supported by experienced technicians and proper equipments, we also provide repairs and calibration for electromedical equipment. ENDO exist to serve your various needs of Medical Equipments and highly concerned about customer satisfaction to promote long term working relationship.



PT. ENDO INDONESIA

Jl. Raya Menganti No. 14, Wiyung - Surabaya 60223

Tel : (031) 767 3636 | Email : info@endo.id

Fax : (031) 767 3737 | URL : <https://endo.id>

Toll Free : 0800 177 ENDO (3636)

