


ENDO

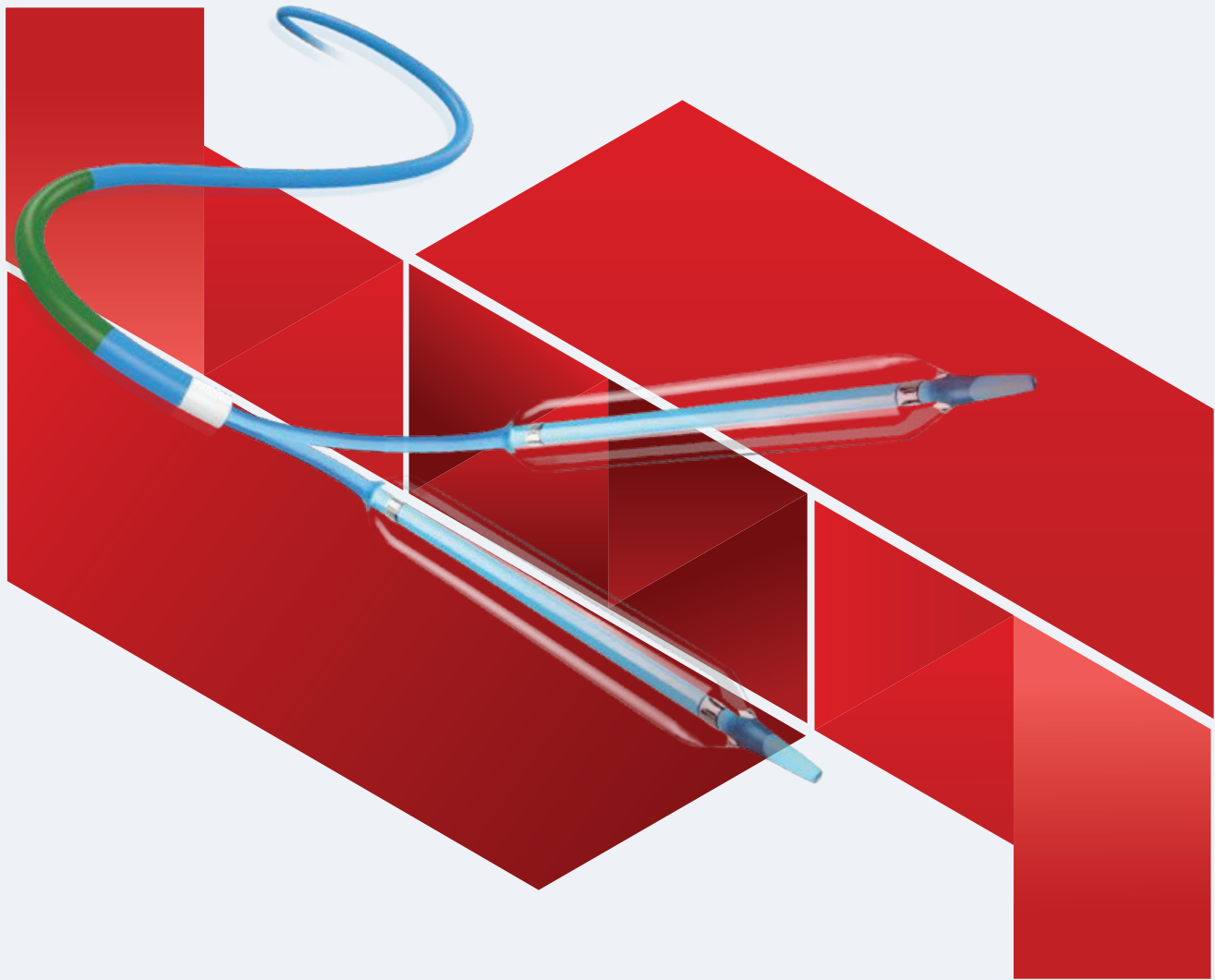

ENDO Mini Coronary Dilatation Catheter

<https://endo.id>

ENDO Mini Coronary Dilatation Catheter

Features:

- Small tip and crossing profiles ensure the better performance for CTO lesions.
- Capable for kissing balloon under 6F guiding catheter.
- Patented hydrophilic coating on distal shaft and Hypotube proximal shaft provide good deliverability.
- Three-directional memory folding of the balloon allows of size minimization with great withdrawal capability, and ensures the good performance of repeated dilatation.
- Soft tapered tip provides smooth transition from guidewire to balloon catheter to avoid “fish mouth”, and minimize the vessel damage.
- Short shoulder design at both ends of the balloon ensures the precise dilatation, and reduces the damage to normal vessel next to the lesion.


Technical Parameters


Catheter Design	Rapid Exchange
Balloon Compliance	Semi-compliant
Marker	2
Distal Shaft O.D.	2.7F
Proximal Shaft O.D.	1.9F
Tip Profile	0.017”
Minimum Guiding Catheter Compatibility	No less than 5F (0.056”)
Guide Wire Compatibility	0.014”
Effective Length	142cm
Nominal Pressure	8 atm (ϕ 2.25mm-3.0mm) / 6 atm (ϕ 3.25mm-4.0mm)
Rated Burst Pressure	16 atm


Ordering Information


Balloon Diameter (mm)	Balloon Diameter (mm)						
	8	10	12	15	20	25	30
1.5	ENDO LPMRX 15008	ENDO LPMRX 15010	ENDO LPMRX 15012	ENDO LPMRX 15015	ENDO LPMRX 15020	ENDO LPMRX 15025	/
1.75	ENDO LPMRX 17008	ENDO LPMRX 17010	ENDO LPMRX 17012	ENDO LPMRX 17015	ENDO LPMRX 17020	ENDO LPMRX 17025	ENDO LPMRX 17030
2	ENDO LPMRX 20008	ENDO LPMRX 20010	ENDO LPMRX 20012	ENDO LPMRX 20015	ENDO LPMRX 20020	ENDO LPMRX 20025	ENDO LPMRX 20030

ENDO

 Jalan Raya Menganti 14 Kedurus, Surabaya 60223
Jawa Timur, Indonesia

 +62-31-767-3636

 info@endo.id

 +62-31-767-3737

 <https://endo.id>

 0800 . 177 . ENDO (3636)