


ENDO

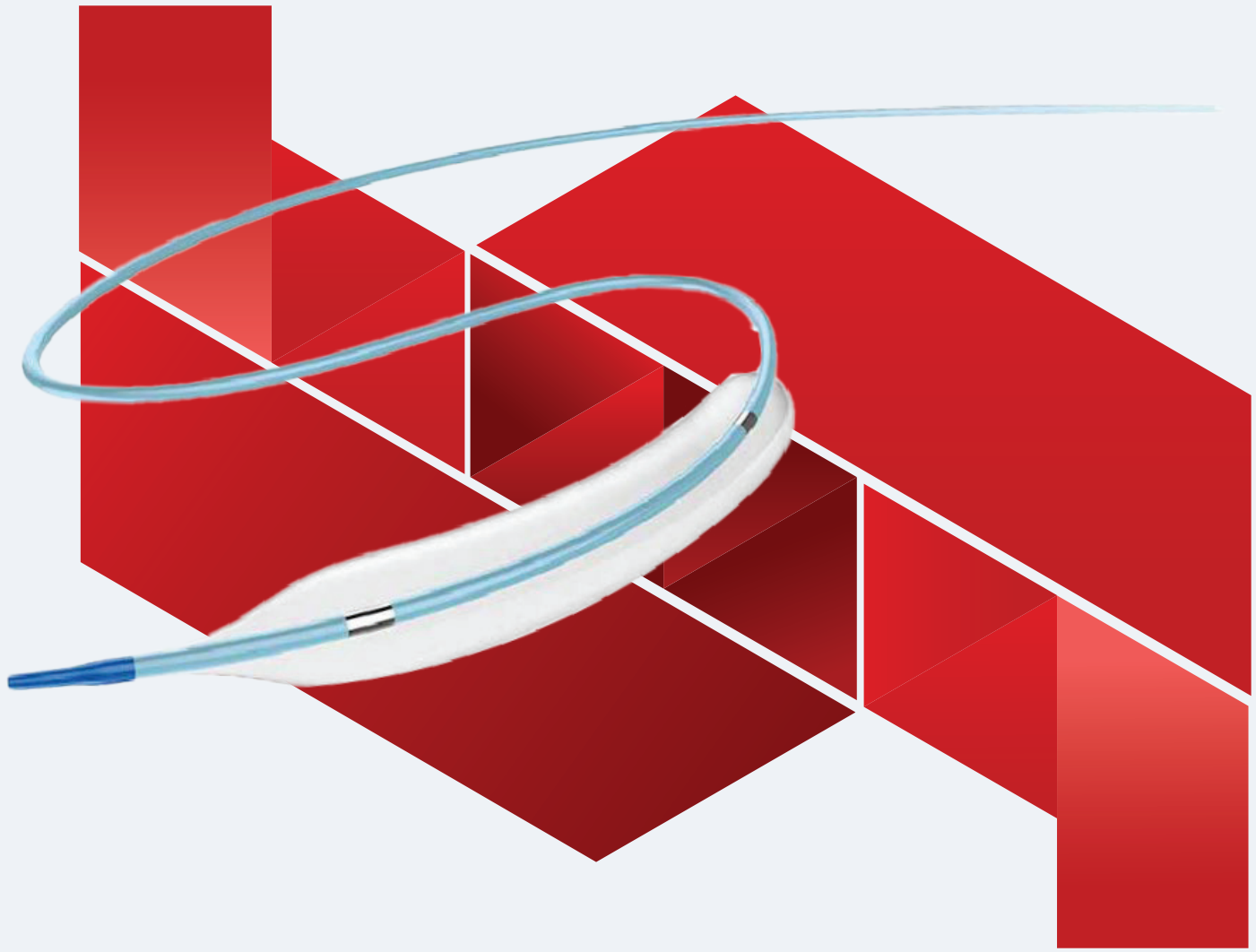

ENDO NC Coronary Dilatation Catheter

<https://endo.id>

ENDO NC Coronary Dilatation Catheter

Features:

- Strong radial strength with high RBP (22 atm) contributes to excellent stent apposition for post-dilatation.
- Patented hydrophilic coating on distal shaft provides good deliverability.
- Three-directional memory folding of the balloon allows of size minimization with great withdrawal capability, and ensures the good performance of repeated dilatation.
- Soft tapered tip provides smooth transition from guidewire to balloon catheter to avoid “fish mouth”, and minimize the vessel damage.
- Short shoulder design at both ends of the balloon ensures the precise dilatation, and reduces the damage to normal vessel next to the lesion.

Technical Parameters


Material	Nylon
Catheter Design	Rapid Exchange
Balloon Compliance	Non-compliant
Marker	2
Distal Shaft O.D.	2.7F
Proximal Shaft O.D.	2.3F
Tip Profile	0.018”
Minimum Guiding Catheter Compatibility	No less than 5F (0.056”)
Guide Wire Compatibility	0.014”
Effective Length	135cm
Nominal Pressure	14 atm
Rated Burst Pressure	22 atm

Ordering Information


Balloon Diameter (mm)	Balloon Length (mm)							
	6	8	10	12	15	20	25	30
1.5	ENDO LPNRX15006	ENDO LPNRX15008	ENDO LPNRX15010	ENDO LPNRX15012	ENDO LPNRX15015	ENDO LPNRX15020	ENDO LPNRX15025	ENDO LPNRX15030
2.0	ENDO LPNRX20006	ENDO LPNRX20008	ENDO LPNRX20010	ENDO LPNRX20012	ENDO LPNRX20015	ENDO LPNRX20020	ENDO LPNRX20025	ENDO LPNRX20030
2.25	ENDO LPNRX22006	ENDO LPNRX22008	ENDO LPNRX22010	ENDO LPNRX22012	ENDO LPNRX22015	ENDO LPNRX22020	ENDO LPNRX22025	ENDO LPNRX22030
2.5	ENDO LPNRX25006	ENDO LPNRX25008	ENDO LPNRX25010	ENDO LPNRX25012	ENDO LPNRX25015	ENDO LPNRX25020	ENDO LPNRX25025	ENDO LPNRX25030
2.75	ENDO LPNRX27006	ENDO LPNRX27008	ENDO LPNRX27010	ENDO LPNRX27012	ENDO LPNRX27015	ENDO LPNRX27020	ENDO LPNRX27025	ENDO LPNRX27030
3.0	ENDO LPNRX30006	ENDO LPNRX30008	ENDO LPNRX30010	ENDO LPNRX30012	ENDO LPNRX30015	ENDO LPNRX30020	ENDO LPNRX30025	ENDO LPNRX30030
3.25	ENDO LPNRX32006	ENDO LPNRX32008	ENDO LPNRX32010	ENDO LPNRX32012	ENDO LPNRX32015	ENDO LPNRX32020	ENDO LPNRX32025	ENDO LPNRX32030
3.5	ENDO LPNRX35006	ENDO LPNRX35008	ENDO LPNRX35010	ENDO LPNRX35012	ENDO LPNRX35015	ENDO LPNRX35020	ENDO LPNRX35025	ENDO LPNRX35030
3.75	ENDO LPNRX37006	ENDO LPNRX37008	ENDO LPNRX37010	ENDO LPNRX37012	ENDO LPNRX37015	ENDO LPNRX37020	ENDO LPNRX37025	ENDO LPNRX37030
4.0	ENDO LPNRX40006	ENDO LPNRX40008	ENDO LPNRX40010	ENDO LPNRX40012	ENDO LPNRX40015	ENDO LPNRX40020	ENDO LPNRX40025	ENDO LPNRX40030
4.5	ENDO LPNRX45006	ENDO LPNRX45008	ENDO LPNRX45010	ENDO LPNRX45012	ENDO LPNRX45015	ENDO LPNRX45020	ENDO LPNRX45025	ENDO LPNRX45030

ENDO

 Jalan Raya Menganti 14 Kedurus, Surabaya 60223
Jawa Timur, Indonesia

 +62-31-767-3636

 info@endo.id

 +62-31-767-3737

 <https://endo.id>

 0800 . 177 . ENDO (3636)