

ENDO Surgical Table Electro Hydraulic ENDO 2000E

Electrohydraulic

Ref: M2C 650 0 0



ENDO Surgical Table Electro Hydraulic ENDO 2000E is designed to accommodate different patient positions required for demanding surgical procedures, with an ultimate focus on patient and user safety. ENDO 2000E Surgical Table is constructed with stainless steel and high quality materials to provide durability and long operating life.

Surgical positions of ENDO 2000E are adjusted by an electrohydraulic system. ENDO 2000E Surgical Table has up/down, sliding, Trendelenburg/reverse, right/left lateral, renal surgery, Trendelenburg and back section upward combination, and back section tabletop positions. Positions are adjusted with remote control, in a quite and nonvibrating manner at constant speed.

ENDO 2000E Surgical Table is equipped with battery recharge system. Movements of the table are fed from battery circuit. Battery power of the table has a 30 minutes operating time in cases of mains failure providing approximately 30 complete position changes.

The tabletop consists of four sections: head, backrest, seat and leg, each with removable, easy to clean, anti-bacterial, and waterproof 80 mm mattress for patient comfort.

Easy to install radiographic tabletop kit* enables using ENDO 2000E for X-ray procedures. Sliding motion of the tabletop enables ENDO 2000E Surgical Table compatible for use with C-Arm radiological equipment.

The durable cast iron table base and electrohydraulic piston is covered with stainless steel. The "T" shape design of the ENDO 2000E allows easy access to the tabletop in any position. Special brakes and anti-static casters allow ease of mobility. Stainless steel side rails accept a variety of available accessories.

ENDO 2000E Surgical Table is a state-of-the-art product and with both quality and competitive pricing it is an excellent choice for today's demanding surgery theaters.

(* *Optional*)

Technical Details


Height adjustment: 725 - 1055 mm
Height adjustment (w/mattress): 805 - 1135 mm
Tabletop: 500 x 1880 mm
Table base: 370 x 830 mm
Head section: 500 x 240 mm
Backrest section: 500 x 490 mm
Seat section: 500 x 510 mm
Right leg section: 235 x 555 mm
Left leg section: 235 x 555 mm
Voltage: 230V, 50Hz
Remote control voltage: DC 5V
Motor inputs: DC 24V, 10A
Battery: 2 each, DC 12V, 7Ah, maintenance free, lead acid battery
Battery recharge voltage: DC 27.2V
Battery recharge time: 18h (full charge)
Low Battery Level Alarm: Sound alarm

Packing Details

Net weight: 250 kg
Gross weight: 295 kg
Packing dimensions: 1000 x 1680 x 1100 mm

Positions

Head section: Upward 15°, Downward 90°, Extendable 5 cm, Removable
Backrest section: Upward 75°, Downward 45°, Upward 45° (w/radiographic tabletop and mattress)
Leg section: Upward 15°, Downward 90°, Spreading 90°, Removable, Independent movement
Lateral: Right 20°, Left 20°
Trendelenburg: 25°
Reverse Trendelenburg: 25°
Lateral nephrectomy: 135°
Reflex: 105°
Flex: 225°
Tabletop: Rotates 355°
Horizontal tabletop sliding: 250 mm

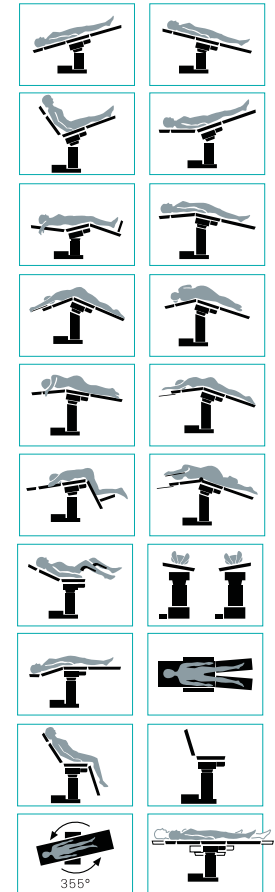

ENDO Surgical Table Hydraulic ENDO 550S

Ref: M2A 000 0 0 - 550S Surgical Table - Hydraulic



ENDO Surgical Table Hydraulic ENDO 550S are designed to accommodate various type of patient positioning required for demanding surgical procedures with ultimate care to the patient's and user's safety. Surgical Tables are constructed with stainless steel and high quality materials to provide durability and long operating life.

The tabletop consists of four HPL-Compact Laminant sections: head, backrest, seat and leg; each with removable, easy to clean, antibacterial, and waterproof polyurethane pads for patient comfort.

The radiographic tabletop* made of HPL - Compact Laminant can be easily installed and removed for X-ray procedures.

The durable cast iron table base is covered with a stainless steel. The "T" shape design of the ENDO 550S allows easy access to the tabletop in any position. Special brakes and antistatic casters allow ease of mobility. Stainless steel side rails accept a variety of available accessories.

ENDO 550S Surgical Tables are reliable products having both quality and competitive pricing it is an excellent choice.

ENDO 550S Surgical Table, Hydraulic

The tabletop can be raised and lowered to the desired height with a hydraulic lift in a quiet and nonvibrating manner by using a foot pedal. Downward speed is adjustable. Other specific table settings are performed manually. The tabletop can be rotated for 355°.

() Optional*

Technical Details


- Height adjustment:** 770 - 1000 mm
- Height adjustment (w/mattress):** 800 - 1030 mm
- Tabletop (w/o accessories):** 450 x 2020 mm
- Table base:** 550 x 800 mm
- Head section:** 450 x 310 mm
- Backrest section:** 450 x 550 mm
- Seat section:** 450 x 445 mm
- Right leg section:** 215 x 680 mm
- Left leg section:** 215 x 680 mm

Packing Details

- Net weight:** 220 kg
- Gross weight:** 265 kg
- Packing dimensions:** 1000 x 1680 x 1100 mm

Positions

- Head section:** Upward 15°, downward 90°, Extendable 5 cm, Removable
- Lateral:** Right 20°, Left 20°
- Trendelenburg:** 25°
- Reverse Trendelenburg:** 25°
- Backrest section:** Downward 15°, Upward 75° (w/mattress), Upward 90° (w/o mattress), Upward 20° (w/radiographic tabletop and mattress)
- Reflex:** 105°
- Flex:** 195°
- Leg section:** Downward 90°, Upward 15° Spreading 90°, Removable, Independent movement
- Renal surgery:** 165°
- Tabletop:** Rotates 355° (550S)

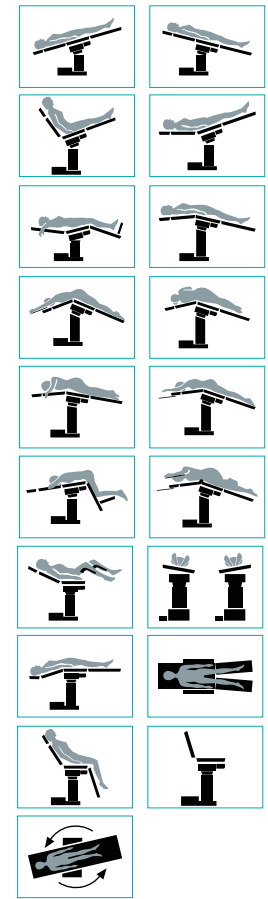

355° (550S)

Accessories For ENDO 2000E, ENDO 550S Surgical Tables

Standard accessories

- ▶ Anaesthetist's frame (1 ea)
- ▶ I.V. stand (1 ea)
- ▶ Shoulder supports (2 ea)
(Optional accessories for 2000E)
- ▶ Body supports (2 ea)
(Optional accessories for 2000E)
- ▶ Surgical table arm boards (2 ea)
- ▶ Knee crutches (2 ea)
(Optional accessory for 2000E,
550S)
- ▶ Compact stand for surgical table
accessories (1 ea)
(Optional accessory for 2000E,
550S)
- ▶ Body strap (1 ea)
- ▶ Mattress (1 set)



Anaesthetist's frame
Ref: M3S 000 0 0



Body supports
Ref: M3R 000 0 0 (Ref of 1 ea)
(Optional accessories for ENDO 2000E)



Body strap
Ref: M9E 000 0 0



I.V. stand
Ref: M3T 000 0 0



Surgical table arm boards
Ref: M3M 000 0 0 (Ref of 1 ea)



Mattress
Ref: M2C 800 5 0 (2000E)



Shoulder supports
Ref: M3P 000 0 0 (Ref of 1 ea)
(Optional accessories for ENDO 2000E)



Knee crutches
Ref: M9N 000 0 0 (Ref of 1 ea)
(Optional accessory for ENDO 2000E,
ENDO 550S)



Polyurethane mattress
Ref: M2C 903 5 0
Ref: M2D 901 5 0 (550)

Accessories

- ▶ Foot support
- ▶ Knee crutches
- ▶ Arthroscopic knee crutch
- ▶ Knee rest
- ▶ Proctology attachment
- ▶ Radiographic tabletop assembly kit
- ▶ Radiographic tabletop kit (550)
- ▶ X-ray cassette holder
- ▶ Wrist grip
- ▶ Arm strap
- ▶ Arm/hand surgical section
- ▶ Waste tray
- ▶ Leg rest - Lithotomy position
- ▶ Head section, high adjustable
- ▶ Arm support
- ▶ Horseshoe headrest
- ▶ Skull clamp
- ▶ Crossbar adaptor for skull clamp
- ▶ Stand for skull clamp and
horseshoe headrest
- ▶ Arthroscopic shoulder support
- ▶ Lithotomy leg positioning boat
- ▶ Spinal surgery frame
- ▶ Carbon fiber leg section (2000E)
- ▶ Compact laminate leg section
(2000E)



Foot support
Ref: M3J 000 0 0 (Ref of 1 ea)



**Compact stand for surgical
table accessories**
Ref: M3E 100 5 0
(Optional accessory for ENDO
2000E, ENDO 550S)



Arthroscopic knee crutch
Ref: M2D 309 5 0



X-ray cassette holder
Ref: M9D 000 0 0



Leg rest - Lithotomy position
Ref: M2D 600 5 0 (Ref of 1 ea)



Horseshoe headrest
Ref: M3F 000 0 0



Arthroscopic shoulder support
Ref: M9S 150 0 0



Knee rest
Ref: M9L 000 0 0 (2000E)
Ref: M30 000 0 0 (550)



Wrist grip
Ref: M2D 333 5 0



Head section, high adjustable
Ref: M9B 200 5 0
Ref: M9B 000 0 0 (550)



Skull clamp
Ref: M9A 000 0 0



Lithotomy leg positioning boat
Ref: AU-OTJ01-001 (Ref of 1 pair)



Proctology attachment
Ref: M3I 000 0 0



Arm/hand surgical section
Ref: M9M 100 0 0



Head section, high adjustable
Ref: M9B 250 0 0 (2000E)



Crossbar adaptor for skull clamp
Ref: M9A 001 5 0



Spinal surgery frame
Ref: M9P 200 0 0



Radiographic tabletop assembly kit
Ref: M2C 400 5 0
Radiographic tabletop kit
Ref: M3V 000 0 (550)



Waste tray
Ref: M9H 010 0 0 (2000E)
Ref: M3Z 000 0 (550)



Arm support
Ref: M9M 000 0 0 (Ref of 1 ea)



Stand for skull clamp and horseshoe headrest
Ref: M9K 000 0 0



Carbon fiber leg section
Ref: M9S 050 0 0 (2000E)



Compact laminate leg section
Ref: M9S 060 0 0 (2000E)