


## PERITONEAL DIALYSIS CATHETER SET CATHETER SET

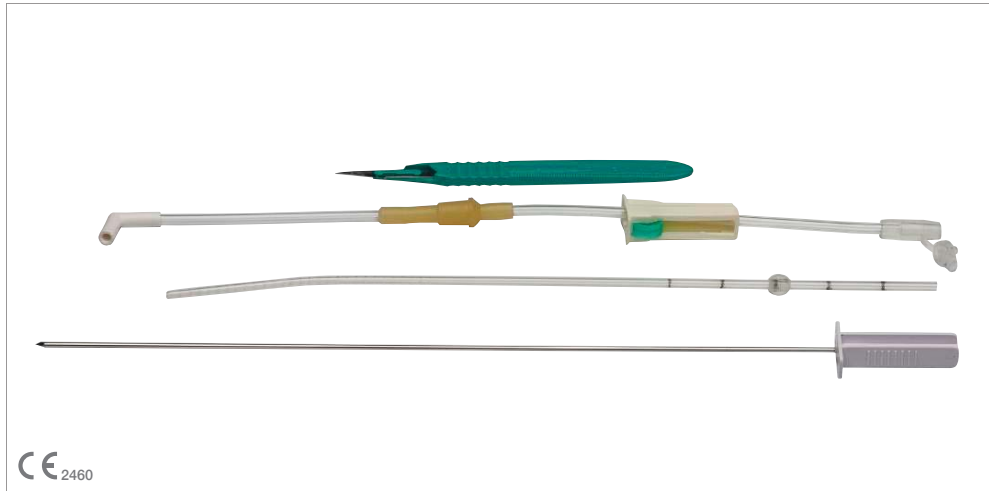

### Peritoneal Dialysis Catheter Set

GS-3010

- Stainless steel trocar is provided with triple angle sharpness to facilitate smooth penetration.
- Spare scalpel blade is supplied to provide the incision.
- Multihole, open end catheter is manufactured from medical grade polyamide.
- The trocar and the tip of the catheter is perfectly matched to facilitate trauma free entry.
- Latex free flash ball type injection port is provided on the junction unit for additional medication.
- Complete kit is placed in thermoformed tray, packed in peelable pouch.
- Sterile, ready for use.
- Box of 20/Master box of 200.

UNIT CONSISTS OF :

- Stainless steel trocar with multihole smooth tip catheter.
- Junction unit with control clamp and scalpel.

### Peritoneal Dialysis Transfusion Set

GS-3012

- Specially designed "Y" shaped transfusion set for administering dialysis solution.
- Provided with two upper control clamps to facilitate alternative change of fluid bottles.
- Two lower control clamps help easy input and drainage of solution.
- Proximal end is provided with suitable adapter for easy connection to catheter unit.
- Sterile, individually packed in peelable pouch.
- Box of 10/ Master box of 100.

#### SPECIFICATIONS

Sizes	Perforation	Marking from Tip
Adult	100 mm	190, 210, 230, 250, 270 mm
Child	50 mm	130, 150, 170 mm

## UMBILICAL CATHETER UMBILICAL VENOUS CATHETER



### Umbilical Catheter GS-3002

- Specially designed catheter for intermittent or continuous access to the umbilical vessels of new born or pre-mature infants.
- Suitable for monitoring blood pressure, withdrawal of blood for sampling without repeated needle pricks, intra venous infusion and exchange transfusion etc.
- Smooth round tip helps easy cannulation.
- Open distal end without lateral eyes, eliminates the chance of clot formation in the blind spaces.
- Radio opaque catheter is provided for X-ray visualisation.
- Clear scale markings at every cm helps to ascertain depth of insertion.
- Proximal end is fitted with female luer mount for convenient connection to I.V therapy equipment.
- Individually packed in peelable pack.
- Sterile, ready for use.
- Box of 20/ Master box of 500.

#### SPECIFICATIONS

Sizes	Colour	Length
3,5	Clear	35 cm
4	Clear	35 cm
5	Grey	35 cm
6	Green	35 cm
8	Blue	35 cm