OUR MISSION

To make a mark in our chosen area, by being ever abreast of new possibilities and tapping the market opportunities.

To anticipate the needs of our customers and end users translating their requirements into meaningful products.

To create, innovate and produce excellence through a philosophy of uncompromising dedication at every level of function.

To make the best use of resources and technology in our pursuit for continuous improvement in our product portfolio.

To nurture a lasting relationship with our distributors - our partners in growth, through empathic understanding and mutual respect.



ABOUT US

One of free India's earliest, most critical revolutions took place in 1952 in the field of healthcare.

At a time when most patients were vulnerable and many succumbed to infections caused by re-usable medical devices, Romsons stepped forward with a pioneering breakthrough - the concept of disposable medical and surgical devices.

Today, over half a century later, the culture of disposable medical devices is well established and Romsons has entrenched itself as the prominent brand in the business -a pioneer and a leader.

An INR 11.0 billion, professionally-managed enterprise.

A global player with presence in over 75 countries.

A retail footprint across geographies - with 1000 plus distributors in India and abroad. Most significantly, the name Romsons continues to inspire trust for unsurpassed quality, innovation and safety of its products.

ROMSONS GROUP OF COMPANIES:

- Romsons Group Private Limited.
- Romsons International-100% Export Oriented Unit
- Romsons Medicons
- Raj Vijay Corporation (Unit I & II)
- RGI Meditech Pvt. Ltd

A PRODUCT FOR EVERY NEED? WELL ALMOST

We have the largest product portfolio in the industry (120 products in over 450 sizes), which caters to almost the entire spectrum of a patient's need. Each product has been designed and manufactured to deliver maximum performance and safety. We have earned recognition in the form of award from The Chamber of Commerce in India for outstanding export performance and the special recognition instituted by Ministry of MSME, Government of India, our pioneering work has been universally appreciated.

WINNING THE TRUST OF THE BEST:

Our experience (over 66 years) in understanding patients' needs and doctors' requirements has earned us the loyalty of the most demanding customers in the world.

Our ability to serve the emerging needs of our customers with speed and agility makes us the preferred choice in the competive market place. Little surprise then, that Romsons products are being exported to developed and emerging markets in over 75 countries.

We also regularly participate in all the major international medical trade fairs like Medica (Germany), Hospitalar (Brazil), Arab Health (UAE), FIME (USA) amongst others.

RESEARCH AND DEVELOPMENT:

A relentless quest: Our mission is to create products that promote healing, well-being and safety of patients. So, we are continuously improving existing products and researching new solutions in the light of changing customer needs and market dynamics. Our research, design and development teams work in tandem to realize the above objective.





GOOD MANUFACTURING PRACTICES:

Our manufacturing plant conforms to GMP, a WHO proprietary benchmark. The discipline of this standard is evident in the physical ambience of our plant, the air and water quality, manufacturing technology, packaging and the sterilization process.

Clean room manufacturing facility ensures minimal microbial/ particle contamination. Employees undergo periodic medical checkup to maintain high standards of health and hygiene. All clean rooms are provided with a three step changing room and air showers to ensure a complete dust-free environment.

MANUFACTURING PLANT:

We have taken great care to source up-to-the-minute, frontline technology. As a result we are able to process moulding and extrusion with precision and within close tolerance.

PACKAGING AND STERILIZATION:

Our products are packed on most modern packaging machines suitable for working in clean room conditions, to ensure minimal risk of contamination or damage during transportation.

We use Ethylene Oxide to sterilize the products, as per ISO 11135 standard.

The entire process is automated and computerized. Each batch is well documented for the sterilization cycle and released for dispatch after written approval of a senior Quality Assurance official.

TOWARDS ZERO DEFECT QUALITY:

Our goal is to manufacture zero-defect products. We have meticulously planned and built a series of multiple online and manual quality control measures into each stage of raw material procurement, manufacturing and packing cycle. Our Quality Control standards conform to ISO 9001: 2015, ISO 13485: 2016 and the European Medical Device Directive 93/42/EEC certifications. Our establishment is also registered with USFDA and MDR 2017/745.



EMPOWERING PEOPLE TO EXCEL:

We provide a professionally nurturant work environment that empowers employees to be self motivated, think out of the box, value excellence and bond with each other to perform as a team.

All across our corporation, our teams function with flexibility, speed and decisiveness in the most challenging environment. We give tremendous freedom and opportunities to all our employees for superior work quality in competitive and ever changing and evolving business environment.

CONTENTS

ANESTHESIA —	
Suction Catheter	Page 10
Guedel Airways	Page 11
Endotracheal Tube	Page 12
Laryngoscope Set	Page 13
Resq Bag	Page 14
Flexi Mask / Flexi Comfort	Page 15
Oxy Set	Page 16
Hi Mask / Oxy Lock	Page 17
Aero Mist / Aero comfort	Page 18
H.M.E Filter / Catheter Mount	Page 19
Epi Kit / Epi Cath	Page 20
SURGERY —	
Romo Vac Set	Page 22
Romo Seal	Page 23
Chest Drainage Catheter	Page 24
Thoracic Trocar Catheter	Page 25
Yankur Suction Set	Page 26
Vaccu Suck Suction Set	Page 27
Yankur Suction	Page 28
Suction Connecting Tubes	Page 30
UROLOGY —	
Urine Collecting Bag	Page 32
Male Cath / Sil Cath	Page 33
Urometer	Page 34
Urometer 500	Page 35
Nelaton Catheter	Page 36
Uro Cath / Silko Cath	Page 37
Mid Stream Urine Collection Set	Page 39

INFUSION —

Intra Cath-2	Page 40
Vennula / Intracan	Page 41
Neo Cath / Neocare / Micron	Page 42
Inflow Set	Page 43
Intra Flow AS	Page 44
Rotaflo	Page 45
Volumetric Plus / Volufix Plus	Page 46
Bi-Valve / Exxacta	Page 47
Vein-O-Line / Vein-O-Line Marvel	Page 48
PM-O-Line / Hi Line	Page 49
Exteena Family	Page 50
Safsite / Hep Cap	Page 51
Centro	Page 52
GASTROENTEROLOGY —	
Ryle's Tube / Levin's Tube / Romolene	Page 54
Salem Sump Tube	Page 55
Infant Feeding Tube	Page 56
MISCELLANEOUS —	
Respirometer	Paae 58
Breazer	_
Peritoneal Dialysis Catheter Set	_
E.C.G Electrodes	
Trimmer Plus	Page 62
Klik Clamp / Cath Clamp	Page 63
Dignity Premium Adult Diapers / Pull-Ups	Page 64
Diginty Drytex / Mattey Underpads	Page 65
Spongee / Sparklean / Kleanzo Wet Wipes	Page 66
Spongee / Poochie Baby Wipes	Page 67
Dispo Guard / Comfit 3D Face Mask	Page 68
Kare Mask / Kare Mask Plus	Page 69
C.V.P Manometer	Page 70



EC CERTIFICATE Full Quality Assurance System

Certificate No.: 11322-2017-CE-IND-NA-PS Rev. 3.0

Project No.: PRJC-195663-2009-PRC-IND Valid Until
22 January 2024

This is to certify that the quality system of:

Romsons International (Unit - II)

59 J(C) & 9, Noida Special Economic Zone, Noida Dadri Road, Phase II, Noida- 201 305, Uttar Pradesh, India

For design, production and final product inspection/testing of: Sterile Disposable Medical Devices

Has been assessed with respect to:

The conformity assessment procedure described in Annex II excluding section 4 of Council Directive 93/42/EEC on Medical Devices, as amended

and found to comply

Further details of the product(s) and conditions for certification are given overleaf.


Place and date: Høvik, 18 May 2021

Check Validity

For the issuing office: Notified Body 2460 DNV Product Assurance AS



Hazem Tinawi Technical Reviewer

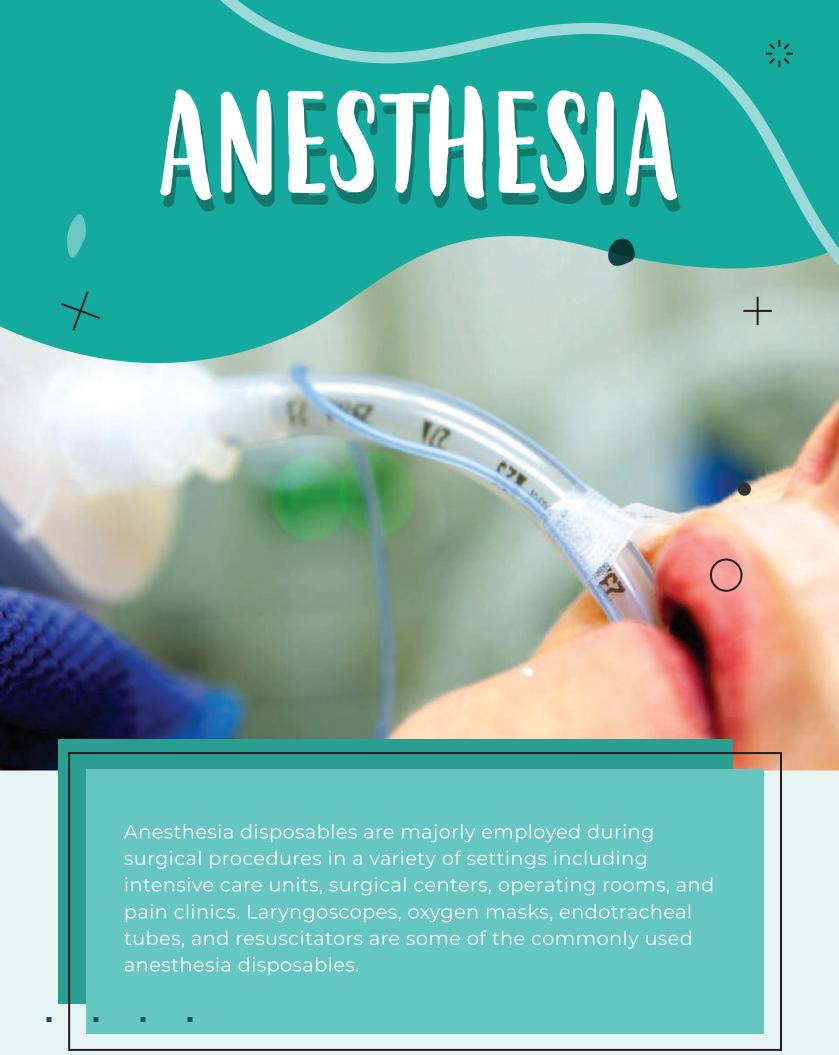

SUCTION CATHETER ENDOBROCHIAL CATHETER



SUCTION CATHETER® GS-2006

- Endo-bronchial catheter to remove respiratory secretions from a patient's airway to keep it clear.
- Suitable for removal of secretions from mouth, trachea and bronchial tubes.
- Provided with universal funnel shaped connector.
- Atraumatic, soft-rounded distal end is open with one lateral eye to avoid trauma during suctioning.
- Frozen surface catheter tubing for smooth intubation.
- · Total length: 53 cm; Effective Length: 50 cm.
- Sterile; Box of 100/ Master Box of 1000.

SUCTION CATHETER® TS WITH STOPPER GS-2030 TS

- Suction catheter with thumb control valve with stopper.
- · Distal open tip with two lateral eyes.
- · Size available from 6 FG onwards.

SUCTION CATHETER® TC GS-2030

 Thumb Control Suction Catheter provided with "T" type vacuum control valve for proper manoeuvring.

SUCTION CATHETER® FT GS-2022

- Finger tip control suction catheter provided with "Y" shape vacuum control valve.
- · Size available from 6 FG onwards.
- · Total length: 53 cm; Effective Length: 47 cm.

SPECIFICATIONS	
Size	Colour Code
5 FG	Grey
6 FG	Sea Green
8 FG	Blue
10 FG	Black
12 FG	White
14 FG	Green
16 FG	Orange
18 FG	Red

GUEDEL AIRWAYS OROPHARYNGEAL AIRWAYS



GUEDEL AIRWAYS GS-2012

- Designed to maintain an unobstructed oropharyngeal airway or following general anaesthesia and in patients who are unconscious for other reasons.
- Manufactured from non-toxic, non-irritant PE polymer.
- Atraumatic soft-rounded edges.
- Provides access for the insertion of suction catheters.
- · Flexible to ensure dental protection.
- Integral hard bite block avoids airway occlusion and biting of tongue.
- Design of the tube conforms to the oropharyngeal anatomy.
- · Stepless airway path for easy cleaning.
- Bite blocks are color coded for instant size identification.
- · Individually packed in peelable soft blister pack.
- Sterile, ready for use.
- Box of 50/ Master Box of 400.

SPECIFICATIONS		
Size	Length	Colour
000	4 cm	Natural
00	5 cm	Blue
0	6 cm	Black
1	7 cm	White
2	9 cm	Yellow
3	10 cm	Red
4	11 cm	Orange
5	12 cm	Purple

ENDOTRACHEAL TUBE PLAIN & CUFFED



ENDOTRACHEAL TUBE PLAIN GS-2002

- Made from thermo-sensitive material readily conforming to body contours at body temperature.
- Kink resistant and transparent material to reduce any possibility of occlusion in the tube bore.
- · Naso-oral tip suitable for nasal/ oral intubation.
- Murphy eye with smooth edges and tip to reduce any trauma to the patient during procedure.
- Radio-opaque for X-ray detection and assistance in placement.
- Proximal end is fitted with standard 15 mm connector which snugly fits with other commonly used connectors.
- · Sterile ; Individually packed in peelable blister pack.
- Box of 10/ Master Box of 200 (size : 2.0 4.5 mm).
- Box of 10/ Master Box of 100 (size : 5 mm).

ENDOTRACHEAL TUBE CUFFED GS-2004

- Thermo-sensitive tube conforms to the tracheal anatomy.
- · Conforms to F-29 standards.
- · Kink resistant tube ensures patient safety.
- Murphy eye has a polished smoothness ensuring atraumatic intubation and extubation.
- · Features a high-volume, low pressure cuff.
- Soft cuff allows perfect seal between tube and trachea.
- Wide black rings proximal to cuff gives visual indication of tube passing the vocal cords to aid apt placement.
- · Radio-opaque line determines exact location.
- · Standard 15 mm connector.
- · Sterile ; Individually packed in blister pack.
- Box of 10/ Master Box of 100.

SPECIFICATIONS	
Product Code	Size (0.5 mm stage)
GS-2002	2 mm - 5 mm

SPECIFICATIONS	
Product Code	Size (0.5 mm stage)
GS-2004	5 mm - 10 mm

LARYNGOSCOPE DISPOSABLE SET



LARYNGOSCOPE SET DISPOSABLE GS-2116

- · Intented for endotracheal intubation.
- · Simple design, textured, sure-grip handle.
- All-in-one, handle and blade combination enables quick size selection from broad range of fibre optic disposable blades.
- Single-Use Blade delivers focused light transmission enabling minimum light scatter and aids clear visualization.
- Disposable use eliminates risk of cross contamination
- · Minimize cost of repeat reprocessing of the device.
- · Reduces intubation time in emergency situations.
- · Type: Macintosh (Curved Blade).
- · Sterile, packed in peelable pack.
- · Master Box of 5/ Master Box of 50.

• SET CONSISTS OF:

- · Laryngoscope Blade
- · Laryngoscope Handle

SPECIFICATIONS		
Fiber Optic Blade	Sizes	

95, 120, 140. 160 mm

Curved

RESQ BAG MANUAL RESUSCITATOR



RESQ BAG® RE-USABLE GS-2050

- Kit consists of a silicone manual resuscitator, star lumen kink resistant oxygen tubing, reservoir bag, pressure limiting valve and silicone mask.
- Made up from 100% silicone.
- Provides bag re-expansion property.
- Pressure limiting valve ensures crucial adjustments for enhanced performance.
- Presence of 360° swivel joint between patient valve and face mask facilitates ease of maneuverability to interface with the patient's posture.
- · Autoclavable manual resuscitator.
- All connectors (15/22mm) as per international standards.
- Box of 1/ Master Box of 8.

SILKO MASK® GS-2049

- Autoclavable mask.
- · Made from 100% Silicone material.
- · Individually packed; Non-Sterile.
- Box of 1/ Master Box of 50.

RESQ BAG® DISPOSABLE GS-2052

- Kit consists of a TPE Manual Resuscitator,
 Multichannel Kink Resistant Oxygen Tubing,
 Reservoir Bag and Inflatable Mask.
- · Resuscitator supports bag re-expansion property.
- Pressure limiting valve ensures crucial adjustments for enhanced performance.
- Presence of 360° swivel joint between patient valve and face mask and facilitates ease of maneuverability to interface with the patient's posture.
- Box of 1/Master Box of 10.

Size
Adult, Child, Infant
0, 1, 2, 3, 4,5

FLEXI MASK/ FLEXI COMFORT OXYGEN MASK



FLEXI MASK® GS-2020

- Designed for oxygen therapy to the patients.
- Moulded face mask manufactured from non-toxic and non-irritant medical grade PVC material.
- · Soft mask fits comfortably over patient's face.
- Adjustable elastic strap for proper positioning of mask on the nasal area.
- 210 cm multichannel, kink-resistant tube ensures continuous oxygen flow.
- Swivel connector for easy attachment to oxygen tube.
- Proximal end fitted with soft funnel shape connector for easy connection to oxygen source.
- · Non-sterile, packed in polybag.
- Individually packed; Master Box of 50.

FLEXI COMFORT TM GS-2401

- Mask with soft seal & strap for better fit & comfort.
- · Manufactured from medical grade PP material.
- · No unpleasant odor.
- Enables the patient to wear the mask for longer duration.
- The cushion also provides better concentration of oxygen as it forms a seal between facial muscles and mask thereby reducing wastage of oxygen.
- Leaves little air blind angle, limiting breath of exhaust.
- Sealing edge of contour line greatly eliminates discomfort caused by the oppression to facial nerve.
- · No metal nose clip, MRI compatible.
- Non-sterile, packed in polybag.
- · Master Box of 50.

SPECIFICATIONS Size Tube Length Adult, Child 210 cm

SPECIAL FEATURES

- ★ Odorless Material ★ No Facial Blind Angle
- ★ Soft Seal at the Edges ★ No Oxygen Wastage

OXY SET TWIN BORE NASAL SET



OXY SET® GS-2016

- Oxy Set is suitable for easy application and efficient administration of oxygen concentration.
- The system of attachment provides maximum freedom to the patient and leaves the patient's mouth free for nutritional and communication purposes.
- Twin prong/ nasal tips are designed to ensure equal volume of oxygen to both the air passages, of the patient.
- Soft funnel shaped connector facilitates easy connection to oxygen source.
- Multichannel tube ensures continuous supply of oxygen, even if the tube is accidentally kinked.
- Manufactured from soft, light weight non toxic material.

AVAILABLE OPTIONS:

- Sterile: Peelable Pack; Box of 20/Master box of 200.
- Non Sterile: Packed in polybag; Master box of 100.

SPECIFICATIONS	
Prong Size	Tube Length
Adult	210 cm
Child	210 cm
Infant	210 cm
Neonate	210 cm

HI MASK/ OXY LOCK OXYGEN MASK



HI MASK® GS-2044

- · Intended for High Concentration Oxygen Therapy.
- Equipped with soft, transparent and odour free mask for greater patient comfort.
- Non-return valves provided at inhale and exhale points, prevents carbon dioxide retention.
- Light weight reservoir improves the oxygen efficiency and patient comfort.
- · Non-sterile, packed in polybag; Master box of 50.
- UNIT CONSISTS OF:
- · Mask with Elastic Strap and Non Return Valve.
- Mask Connector with Non Return Valve.
- · 600 ml Reservoir capacity.
- 210 cm. (7 ft.) long Multichannel Tube.

OXY LOCK® GS-2038

- Variable venturi system to ensure accurate concentrations of oxygen with single unit.
- Transparent air entrainment mask with soft anatomical form with easily adjustable elastic strap.
- Wide and adjustable strap for proper positioning of mask on the nasal area.
- · Rotating type connector for patient's comfort.
- Unique spring loaded mechanism locks prescribed concentration in place ensuring safety & accuracy.
- 210 cm (7 ft.) multichannel tube provided to ensure continuous flow of oxygen.
- Non-sterile, packed in polybag; Master box of 50.

SPECIFICATIONS	
Reservoir	Tube Length
600 ml	210 cm

SPECIFICATIONS		
Mask Size	Tube Length	
Adult, Child	210 cm	

AERO MIST/ AERO COMFORT AEROSOL SET



AERO MIST® GS-2074

- · Soft, clear aerosol mask with anatomical form.
- Gently rolled, feathered edges with nose bridges for extra comfort.
- · Rotating type connector for patient's comfort.
- The kit has nebuliser chamber to house the drug in liquid form with a funnel shaped vent to allow high flow of oxygen converting the liquid drug into fine mist to reach lungs.
- Larger surface area provided by unique convex cone design, ensures maximum capillary action and eliminates medication wastage.
- Nebulizes 2 cc within 10 minutes in horizontal or vertical position ensuring patient comfort.
- Non-sterile, packed in polybag; Master box of 50.

AERO COMFORT TM GS-2403

- · Mask with soft seal & strap for better fit & comfort.
- · Manufactured from medical grade PP material.
- · No unpleasant odor
- · Enables the patient to wear the mask for longer.
- The cushion also provides better concentration as it forms a seal between facial muscles and mask thereby reducing the wastage of medication.
- Leaves little air blind angle that limits breath of exhaust.
- Sealing edge of contour line greatly eliminates discomfort caused by oppression to facial nerve.
- No metal nose clip, MRI compatible.
- · Available Sizes : Adult and Child,
- Tube Length: 210 cm.
- Non-sterile, packed in polybag; Master box of 50.

SPECIFICATIONS Size Tube Length Adult, Child 210 cm

SPECIAL FEATURES

- ★ Adjustable Jet Nebulizer ★ No Medication Loss
- ★ 3cc Nebulization in 10mins ★ Fine Aerosol Mist

H.M.E FILTER/ CATHETER MOUNT ACCESSORIES



H.M.E FILTER GS-2095

- Facilitate to retain heat and moisture from the patient's exhaled gas.
- High efficiency electrostatic filters provide effective protection against various types of particles including bacteria, viruses, and moisture droplets.
- Protect the patient, the equipment, and the breathing circuit from any kind of contamination.
- Gives a low and stable breathing resistance even when in wet condition.
- · Conveniently located port permits gas sampling.
- Standard connections compatible with different breathing systems.
- · Sterile, individually packed.
- Box of 10/ Master Box of 200.

CATHETER MOUNT GS-2404/2407

- Used as an intermediary connection between the patient and the breathing system.
- · Corrugated tube made up of medical grade PE.
- Double flip top cap with seal having provision for both bronchoscopes and for suction catheters.
- Double swivel elbow with 15 mm I.D. / 22 mm
 O.D. connector fits all standard.
- · tracheal tubes & all ventilator circuits.
- · Sterile, individually packed.
- Box of 10, Master box of 200.
- Available in two versions :
- · GS-2404 : Flip Tap[™] Catheter Mount[®]
- CS-2407 : Fip Tap[™] Tip Up Catheter Mount[®]
 (Collapsible)

SPECIFICATIONS		
Size	: Tidal Volume	: Moisture Output
Adult	150-1500 ml	33.5 mg H₂O/ L
Child	70-250 ml	32.3 mg H₂O/ L
Crind	. 70 230 1111	. 32.3 mg 11 ₂ 0/ E

SPECIFICATIONS	
Port Type	Diameter
Bronchoscopy	3.40 mm
Suction	9.40 mm

EPIDURAL KIT EPIDURAL ANESTHESIA KIT




EPI KIT® GS-2026

- KIT CONSISTS OF:
- EPIDURAL/TUOHY NEEDLE: Tuohy needle provided with clear depth marking for accurate insertion and polished inner bevel edge of the needle to minimize the risk of catheter shearing.
- EPIDURAL CATHETER: Closed round catheter tip with three lateral eyes.
- LOR SYRINGE: Specially designed 'loss of resistance' syringe suitable for use with air or saline technique; provides clear identification of epidural space.
- THREADING ASSIST GUIDE: Fits into the edpidural needle and assists in easy advancement of catheter.
- DETACHABLE CONNECTOR: Separate twist and lock connector for connection.
- FLUID FILTER: 0.2 micron filter provides additional safety and control against microbial infections.
- FILTER FIXATOR: To secure fluid filter in place.
- TRANSPARENT DRESSING: To secure the catheter.
- Box of 1/Master Box of 50.

EPI CATH[®] GS-2018

- Epidural catheter is specially designed for continuous anaesthesia and pain relief therapy.
- Catheter is made from bio-compatible plasticizer free polymer.
- Flexible tip with three lateral eyes.
- Catheter is marked to determine the depth of insertion and accurate placement of catheter in epidural space.
- · Luer lock twist detachable connector.
- · Crystal clear catheter with radio opaque line.
- · Sterile, individually packed in peelable pack.
- Box of 10/Master box of 400.

SPECIFICATIONS	
Product Code	Catheter Size
GS-2026/2018	Passes through NG 16,18

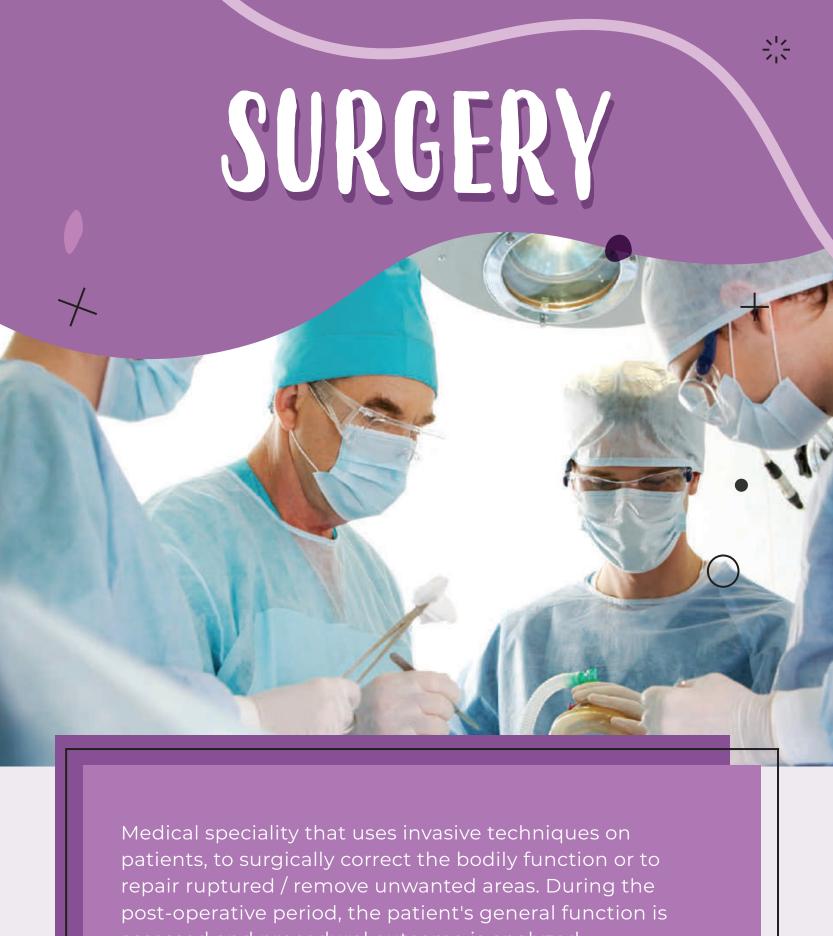

assessed and procedural outcome is analyzed.





ROMO VAC SET CLOSED WOUND SUCTION SET



ROMO VAC® SET GS-5002

- Close wound drainage device under negative pressure.
- Provided with the option to use one or two catheters post-operatively.
- Kink resistant connecting tube.
- Multiperforated catheters provided with radio opaque line and satin smooth perforation/ eyes for trauma free performance.
- Available in different sizes to meet moderate to heavy drainage needs.
- Graduated bellow allows user to measure the drained volume.
- Flexible bellow chamber is easy to depress for activating negative ressure.
- · Sterile ; Individually packed.
- Box of 1/ Master Box as per size.

UNIT CONSISTS OF:

- · Bellow capacity 800 ml/500ml/400ml/50ml.
- · Connecting tube with clamp and "Y" connector.
- · Redon (curved needle) with matching catheter.
- · Spare perforated catheter.

SPECIFICATIONS	
Size	Colour Code
8 FG 10 FG 12 FG 14 FG 16 FG 18 FG 20 FG	Blue Black White Green Orange Red Yellow
Bellow Capacity	Master Packing
GS-5002 S : 800 ml GS-5002 F : 500 ml GS-5002 M : 400 ml GS-5002 N : 50 ml	20 40 40 50

ROMO SEAL UNDERWATER SEAL DRAINAGE SYSTEM




ROMO SEAL[®] GS-5020/ 5031/ 5032

- · Under water seal drainage system.
- Specially designed drainage bottle suitable for pleural drainage in conjunction with intercoastal drainage catheters.
- Double chamber and single chamber compact unit suitable for all cardiac and thoracic procedures.
- Range suitable for management of major, moderate and minor drainages.
- Easy to read markings help note the drained volume.
- Clearly marked initial level ensures effective under water seal.
- Write on facility on the side of scale helps in monitoring the collected drained volume precisely.
- Suction port is provided for connection with suction unit.
- Kink resistant large bore tubing facilitates unrestricted flow.
- Metallic hangers for easy attachment to patient's bedside.
- Sterile, individually packed.
- Box of 1/ Master Box of 10.

SPECIFICATIONS

Product Code	Size	Capacity
GS-5020	Adult	2000 ml
GS-5031	Midi	1200 ml
GS-5032	Kid	500 ml

CHEST DRAINAGE CATHETER INTERCOSTAL

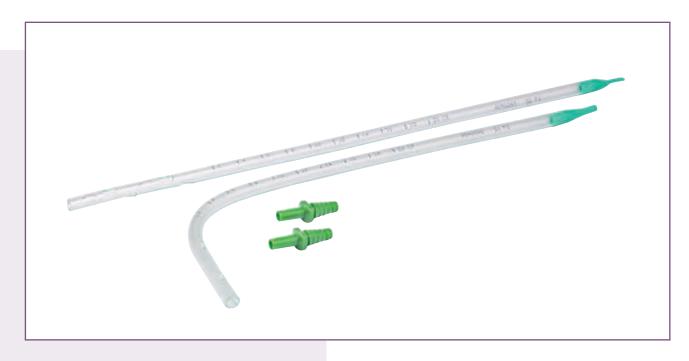

CHEST DRAINAGE CATHETER® GS-5010

- Suitable for intercostal drainage post cardiothoracic surgery.
- Distal end is open with super smooth large eyes to maximize the drainage from the chest.
- · Dark green line for x-ray detection.
- Matching size loose connector facilitates easy connection to drainage system.
- Catheter is marked at every 2 cm to ascertain placement.
- · Sterile ; Box of 20/ Master box of 200.

FLEXO CATH[®] GS-5034

Thoracic drainage catheter, made with extra soft DEHP free material that is gentle to the body tissue.

CHEST DRAINAGE CATHETER® ANGLED GS-5011

- Specially designed right angled catheter for thoracic drainage from the deeper sites of the pleural cavity.
- Dark green line for x-ray detection.
- Matching size loose connector facilitates easy connection to drainage system.
- Sterile, individually double packed in a peelable pouch.
- Box of 10/ Master box of 100.

SPECIFICATIONS		
Size	Length	
12, 16, 20, 24, 28, 32, 36 FG	450 mm	

SPECIFICATIONS	
Size	Length
12, 16, 20, 24, 28, 32, 36 FG	450 mm

THORACIC TROCAR CATHETER TROCAR CATHETER




THORACIC TROCAR CATHETER® GS-5030

- Specially designed for quick non-operative pleural and chest drainage for the relief of chest conditions such as tension pneumothorax.
- Manufactured from medical grade polymer fits snugly over the trocar.
- · Blunt yet smooth tip.
- · Flat handle facilitates easy manoeuvrability.
- · Large, smooth lateral eyes for efficient drainage.
- Markings provided at 5, 10, 15 and 20 cm from the distal tip help ascertain depth of penetration in the intercoastal space.
- Proximal end is fitted with tapered connector for easy connection to the drainage bottle.
- · Individually double packed in peelable pouch.
- Sterile; Box of 10/ Master box of 100.

SPECIFICATIONS	
Size	Length
8 FG	80 mm
10 FG	160 mm
12 FG	160 mm
14 FG	160 mm
16 FG	190 mm
18 FG	350 mm
20 FG	350 mm
24 FG	350 mm
28 FG	350 mm
32 FG	350 mm
36 FG	350 mm

YANKAUR SUCTION SET YANKAUR SUCTION

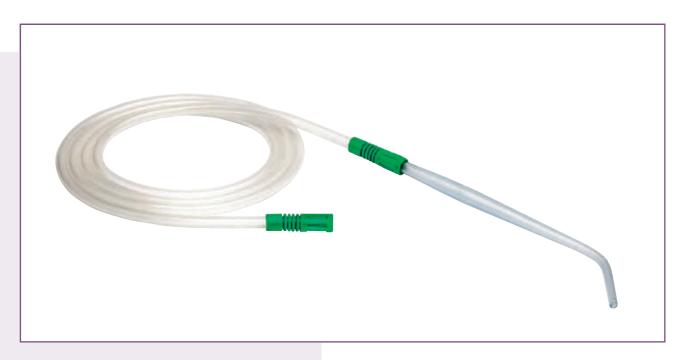

YANKAUR SUCTION SET GS-5016 A6

- Yankaur Suction Set is suitable for convenient removal of secretions, blood & debris etc perioperatively.
- Complete set is provided with Yankaur suction tip mounted on 2 mtr long ribbed tube, provided with universal connector at both ends.
- Universal connectors are moulded from soft PVC, to accommodate the suction tip on one end and all type of connectors of suction apparatus at the other end.
- Kink resistant ribbed tube has the strength to withstand the suction.
- · Individually packed in peelable soft blister pack.
- · Sterile, ready for use.
- Box of 10/ Master box of 100.

SPECIFICATIONS	
Tube Length	I.D/ O.D
200 cm	6/9 mm

VACCU SUCK SUCTION SET YANKAUR SUCTION



VACCU SUCK SUCTION SET® GS-5016

- Moulded Yankaur suction handle attached with 2.5 mtr long tube.
- Flexible, kink resistant tube provides unobstructed suction during the prolonged use.
- Soft flexible adaptors at both ends of the tube provides safe grip to the handle as well as at the suction source.
- Vent port is permanently closed with tight fit sleeve, which can be removed with small incision by blade to change over to vent control system.
- Sterile, Box of 10/Master box of 100.
- Available in two versions :

GS-5016 S: Handle with Standard Tip GS-5016 C: Handle with Crown Tip

SPECIFICATIONS	
Tube Length	I.D/ O.D
250 cm	6/ 9 mm

YANKAUR SUCTION HANDLE GS-5014

- Transparent moulded handle for peri-operative suction.
- Suitable for use with or without vent control facility for better suction management.
- Vent port is permanently closed with tight fit sleeve, which can be removed with small incision by blade to change over to vent control system.
- · Light in weight.
- · Sterile, individually packed.
- Box of 25/Master box of 250.

• Available in two versions :

GS-5014 S: Standard tip with four side holes GS-5014 C: Crown tip with four side holes

YANKAUR SUCTION SET YANKAUR SUCTION

PRODUCT CODE	SUCTION HANDLE / TIP	PRODUCT CODE	SUCTION SETS
GS-5014S	YANKAUR SUCTION HANDLE (moulded) standard open tip With 4 Holes. I.D. / O.D. : 3.75 / 7.30mm	GS-5016S	VACCU SUCK SUCTION SET with standard open tip handle With 250cm Tube 5014S Handle.
GS-5014C	YANKAUR SUCTION HANDLE (moulded) crown tip With 4 Holes. I.D. / O.D. : 3.90 / 10.75mm	GS-5016C	VACCU SUCK SUCTION SET with crown tip handle With 250cm Tube 5014C Handle.
GS-5014 A1 GS-5014 A2	YANKAUR SUCTION TIP - FINE without vacuum control. YANKAUR SUCTION TIP - FINE with vacuum control Length 250mm - With 2 Holes. I.D. / O.D. : 2.50 / 4.00mm	GS-5016 A2	YANKAUR SUCTION SET with fine vacuum control tip With 200cm Long Tube

YANKAUR SUCTION YANKAUR SUCTION

PRODUCT CODE	SUCTION HANDLE / TIP	PRODUCT CODE	SUCTION SETS
GS-5014 A3	YANKAUR SUCTION TIP - MINI without vacuum control. YANKAUR SUCTION TIP - MINI with vacuum control Length 160mm - Without Hole. I.D. / O.D. : 3.20 / 4.80mm	GS-5016 A4	YANKAUR SUCTION SET with mini vacuum control tip With 200cm Long Tube
GS-5014 A5 GS-5014 A6	YANKAUR SUCTION TIP - MIDI without vacuum control. YANKAUR SUCTION TIP - MIDI with vacuum control Length 250mm - With 4 Holes. I.D. / O.D. : 4.00 / 6.00mm	GS-5016 A6	YANKAUR SUCTION SET with midi vacuum control tip With 200cm Long Tube
GS-5014 A7 GS-5014 A8	YANKAUR SUCTION TIP - MAXI without vacuum control. YANKAUR SUCTION TIP - MAXI with vacuum control Length 250mm - With 4 Holes. I.D. / O.D. : 5.20 / 8.00mm	GS-5016 A8	YANKAUR SUCTION SET with maxi vacuum control tip With 200cm Long Tube

SUCTION CONNECTING TUBES CONNECTING TUBE



SUCTION CONNECTING TUBE TM GS-5017


- Suction connecting tube with Moulded and flexible connectors ensure secure connection.
- · Sterile, Individually packed.
- Box of 10/ Master Box of 100.
- Available in :

GS-5019: Suction Connecting Tube Integral Funnel Ends.

CONDUCTIVE SUCTION TUBE ™ GS-5101

- Suction connecting tube with conductive line to prevent static charge.
- Clear, non-coiling tube's wall thickness prevents it from collapsing when under high negative pressure.
- · Sterile, Individually packed in a peelable pouch.
- · Non-Sterile rolls with bubble at each 100 cm.

SPECIFICATIONS	
Tube Length	I.D/ O.D
GS-5017 : 200, 250, 300, 400 cm	Tube : 6 / 9 mm
GS-5019 : 180, 220,	Tube : 7 / 10 mm
270, 300 cm	Ends : 10 / 15 mm
Optional : 30 mtr	Tube : 7 / 10 mm
Non-Sterile Roll	Ends : 10 / 15 mm
GS-5101 : 180, 220,	Tube : 5 / 7 mm
270, 300 cm	Ends : 10 / 13.5 mm
Optional : 30 mtr	Tube : 5 / 7 mm
Non-Sterile Roll	Ends : 10 / 13.5 mm

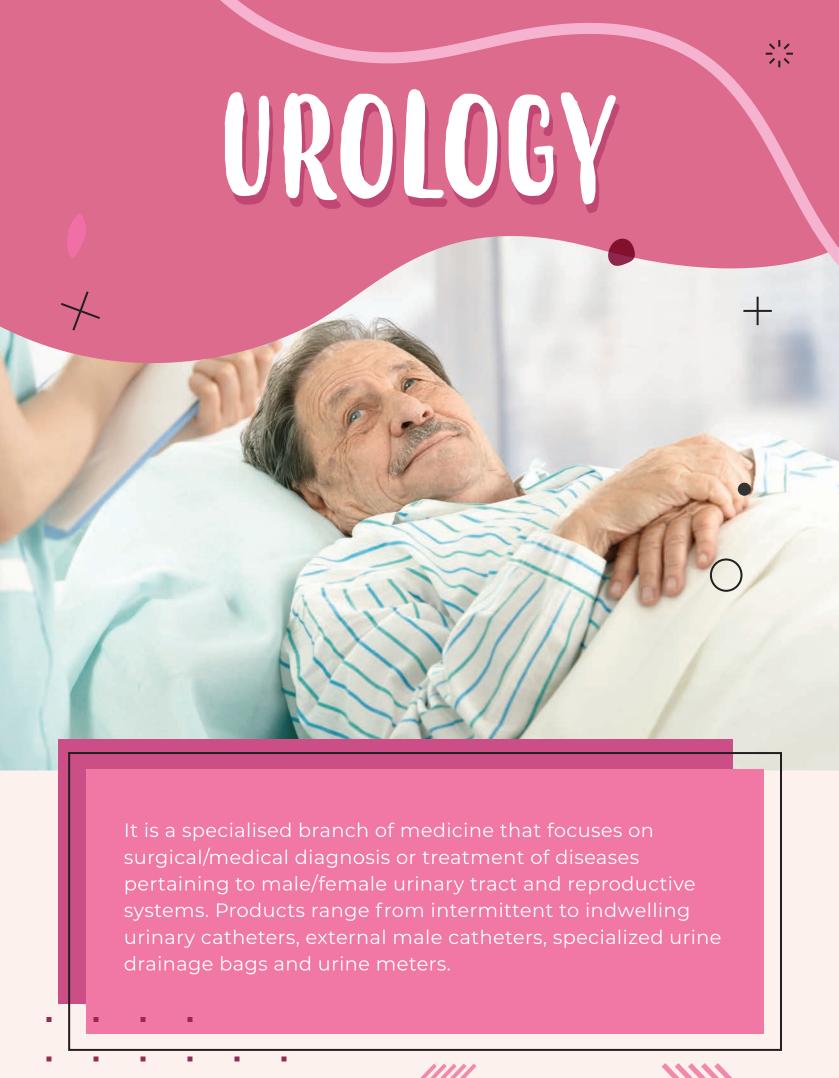

URINE COLLECTING BAG URINE BAG



URINE COLLECTING BAG GS-1020

- · Closed urine drainage system.
- Made from soft clinical grade PVC sheeting.
- Reinforced eyelets allow easy hanging of the bag and provide maximum protection against tearing.
- Push/Pull type bottom outlet allows convenient emptying.
- Kink resistant 90 cm tube provided with universal tapered connector.
- Easy to read scale for quick determination of urine volume.
- Sterile individually packed, Pack of 25 / Master box of 250
- Non sterile bulk pack, Box of 200 / Master box of 400.

AVAILABLE IN FOLLOWING VERSIONS:

- 1020A: Non return valve & bottom outlet.
- 1020AH: Non return valve, bottom outlet & carrying hanger.
- 1020T: Non return valve & 'T' tap bottom outlet.
- 1020B: Non return valve & without outlet.

SPECIFICATIONS	
Capacity	Tube Length
2000 ml	90 cm

MALE CATH/ SIL CATH EXTERNAL CATHETER



MALE CATH® GS-1010

- Male Cath is specially designed for urine incontinence for day/night use in male patients.
- · Made from 100% latex for a soft and gentle feel.
- Provided with soft adhesive coated strip for proper fixation of the catheter.
- Proximal end is tubular shaped, designed for safe and easy connection to the urine bag or leg bag.
- · Individually packed in peelable soft blister pack.
- Sterile, ready for use.
- Box of 50/Master box of 1000.

SPECIFICATIONS Size Tube Length Small 20 mm Medium 25 mm Large 30 mm Extra Large 35 mm

SIL CATH[®] GS-1027

- · 100% silicone male external catheter.
- Soft silicone catheter is designed to be nonirritating, non-reactive for a soft and gentle feel.
- · Cup-shaped applicator makes it easy to wear.
- Ergonomically shaped and easily rolls over the penis.
- Transparent material allows better skin health monitoring while the catheter is in place.
- Provided with soft adhesive strip for proper fixation of the catheter to penis shaft.
- Skin-friendly adhesive helps to maximize wear time and also prevents accidental detachment.
- Gentle adhesive keeps the catheter in place yet allows for easy removal.
- Proximal end is designed for safe and easy connection to urine bag/ leg bag.
- Sterile, individually packed in peelable soft blister pack.
- Box of 25/ Master Box of 400.

UROMETER MEASURED CHAMBER WITH URINE BAG




UROMETER® GS-1018

- · Designed for measurement of hourly urine output.
- Clear and transparent measured volume chamber provides accurate measurement of the urine.
- Urine chamber directly attached to the urine bag with overflow facility.
- 150 cm long kink resistant tubing ensures unrestricted urine flow.
- Completely closed system eliminates the risk of contamination.
- Screw type bottom outlet provided for easy and quick emptying of urine bag.
- Universal tapered connector provided with sampling port facilitates mid-stream urine sampling.
- · Sterile and individually packed in a box.

SPECIFICATIONS		
Product Code	Box/ Master Box	
GS-1018 500E	10/50	
GS-1018 250	5/50	
GS-1018 100	5/50	

SPECIFICATIONS		
Product Code	Chamber Cpty.	Bag Cpty.
GS-1018 500E	500 ml	2000 ml
GS-1018 250	250 ml	2000 ml
GS-1018 100	100 ml	800 ml

UROMETER 500 MEASURED CHAMBER WITH URINE BAG

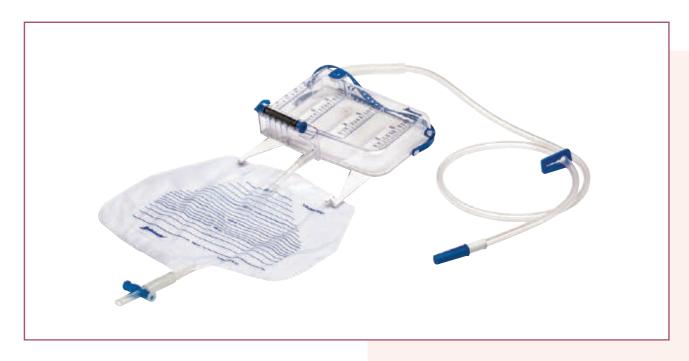


UROMETER[®] **500** GS-1018-500

- Closed urine drainage system with graduated chamber to measure the urine output on hourly basis in ICU and CCU's.
- 500 ml transparent chamber has four fine measuring sections to measure the accurate and precise volume of the urine.
- Collection chamber is provided with hydrophobic filter to minimize the risk of urinary tract infection.
- Single handed operation of valve for emptying all chambers simultaneously.
- Bypass system is provided to prevent the overflow.
- 2 litre capacity urine collecting bag with clamp and "T" shaped bottom outlet.
- 120 cm long kink resistant tubing ensures unrestricted flow of urine.
- Tapered connector provided with sampling port facilitates easy and aseptic mid-stream urine sampling.
- Flexible fastening strap for easy hanging to bed side.
- · Sterile, individually packed in box.
- Box of 1/Master box of 20.

SPECIFICATIONS	
Chamber Capacity	Bag Capacity
500 ml	2000 ml

NELATON CATHETER URINE DRAINAGE CATHETER

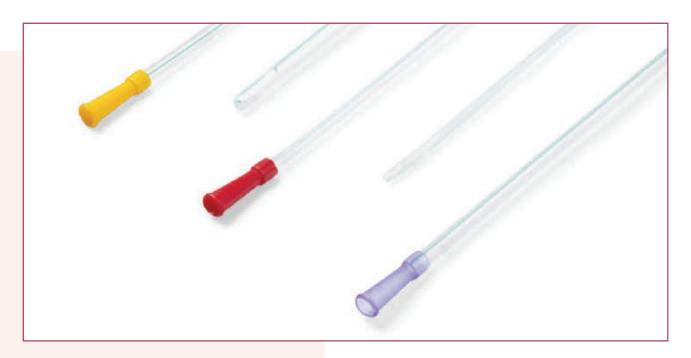

NELATON CATHETER GS-1004

- Nelaton catheters are used for short term bladder catheterisation.
- Manufactured from non-toxic, non-irritant medical grade PVC compatible with catheter lubricants.
- Frozen surface tubing for smooth intubation.
- Presence of radio-opaque line throughout the catheter tube.
- Distal end is coned for non-traumatic introduction with two lateral eyes for efficient drainage.
- Proximal end is fitted with universal funnel shaped connector for extension.
- Funnel shape connectors are colour coded for instant identification of sizes.
- · Straight packed in peelable soft blister pack.
- Sterile; Box of 100/ Master box of 1000.

• AVAILABLE IN FOLLOWING VERSIONS :

GS-1004 : MALE

• Effective Length: 37cm, Total Length: 40 cm.

GS-1004 F : FEMALE

• Effective Length: 17cm, Total Length: 20 cm.

SPECIFICATIONS	
Size	Colour Code
6 FG	Sea Green
8 FG	Blue
10 FG	Black
12 FG	White
14 FG	Green
16 FG	Orange
18 FG	Red
20 FG	Yellow
22 FG	Violet
24 FG	Blue

URO CATH FOLEY BALLOON CATHETER



URO CATH® 2 WAY GS-1072

- · Foley catheter used for urine drainage.
- Made up of natural rubber latex bonded with silicone elastomer which eliminates the risk of encrustation.
- · Smooth surface facilitates non-traumatic intubation.
- Distal end is coned and provided with super smooth, burr free eyes for atraumatic intubation.
- Provided with symmetrical large capacity balloon that ensures straight tip for proper urine flow and good sphincter action, preventing bladder leakage.
- Hard non-return valve ensures easy inflation and deflation of balloon.
- · Valve is colour coded for instant size identification.
- Sterile; Box of 50/ Master Box of 500.

AVAILABLE IN FOLLOWING VERSIONS :

• GS-1094 : URO CATH[®] 3 WAY

Size : 16 - 24 FG

• GS-1078 : SILKO CATH[®] (100% Silicone Catheter)

• Size: 8 - 24 FG; Box of 10/ Master Box of 500.

SPECIFICATIONS

Size	Balloon Capacity
6 FG	1.5- 3 ml
8 FG	5 ml
10 FG	5 ml
12 FG	5-15 ml
14 FG	30 ml
16 FG	30-50 ml
18 FG	30-50 ml
20 FG	30-50 ml
22 FG	30-50 ml
24 FG	30-50 ml

MIDSTREAM URINE COLLECTION COLLECTION SET




MIDSTREAM URINE COLLECTION SET (WITH FUNNEL) GS-252

- Device Urine Collection Set intended for midstream sampling.
- Consists of a transparent conical based container with easy to use funnel and a screw fit cap.
- Funnel facilitates easy collection of midstream urine samples.
- Device helps reduce any risk of microbial contamination.
- Write-on facility on the transparent container to mention patient details.
- Container's conical base encourages quick visualization of urine sediments.
- · Sterile, individually packed.
- · Container capacity: 20 ml
- Master Box of 100

• AVAILABLE IN FOLLOWING VERSIONS:

 GS-257: Mid Stream Urine Collection Set (with Funnel and Non-Woven Wipe)

SPECIFICATIONS	
Product Code	Capacity
GS-252	20 ml
GS-257	20 ml

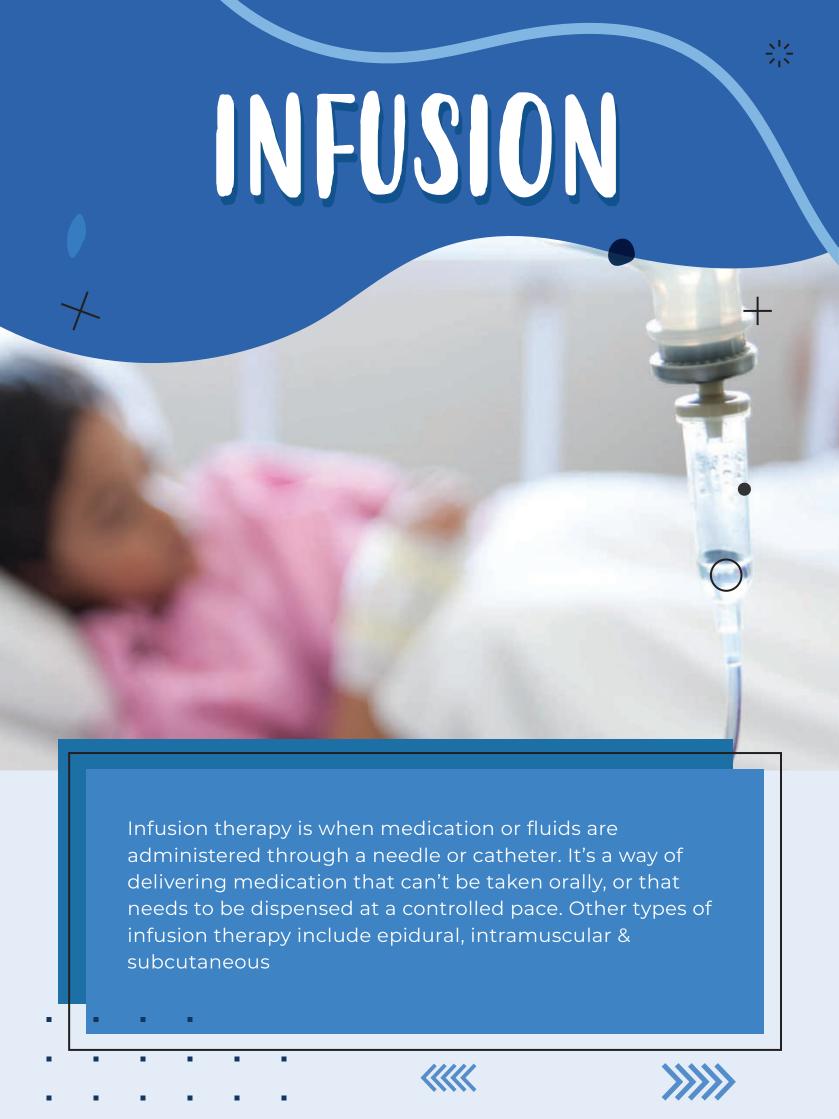

INTRA CATH-2 INTRAVENOUS CANNULA




INTRA CATH®-2 GS-3022

- I.V Cannula with wing, injection port and luer lock.
- Manufactured from fully automated CAD/ CAM process.
- Super sharp triple facetted beveled needle is designed with two cutting edges for smooth penetration.
- · Kink resistant FEP-Teflon cannula
- Extremely smooth inner and outer surface of cannula ensures smooth flow.
- · Double tapered molded cannula tip.
- 3 Radio opaque stripes cannula for X-ray visualization.
- · Transparent flashback chamber for blood detection.
- · Individually packed in blisted pack.
- Sterile; Box of 100/ Master Box of 1000.

SPECIFICATIONS			
Size	Length	Colour	
14 G	45 mm	Orange	
16 G	45 mm	Grey	
18 G	45 mm	Green	
20 G	32 mm	Pink	
22 G	25 mm	Blue	
24 G	19 mm	Yellow	

VENNULA/ INTRACAN INTRAVENOUS CANNULA

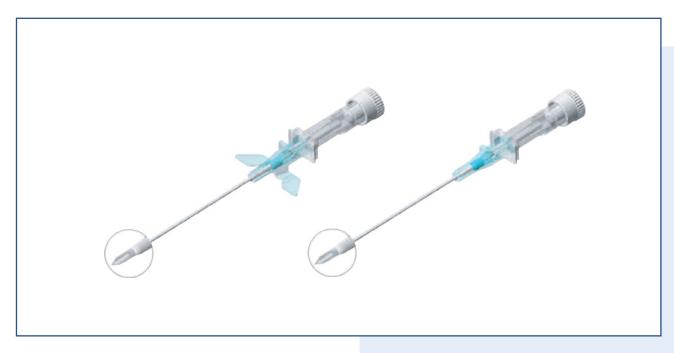

VENNULA® GS-3033

- I.V Cannula without injection port and without wings.
- Kink resistant FEP-Teflon cannula has extremely smooth inner/ outer surface ensures minimum trauma and smooth flow.
- Super sharp triple facette beveled needle is specially designed with two cutting edges for smooth penetration.
- Transparent flashback chamber allows immediate detection of blood.
- · Available With Filter & Without Filter.
- Sterile; Box of 100 / Master box of 1000.

INTRACAN® GS-3225

- · I.V Cannula with wings and without injection port.
- Small, flexible wings and absence of injection port permits easy handling, maneuvering and secure fixation.
- · Transparent flashback chamber.
- Available Sizes: 18, 20, 22, 24 FG

SPECIFICATIONS			
Size	Length	Colour	
14 G	45 mm	Orange	
16 G	45 mm	Grey	
18 G	45 mm	Green	
20 G	32 mm	Pink	
22 G	25 mm	Blue	
24 G	19 mm	Yellow	

NEO CATH/ MICRON PAEDIATRIC I.V CANNULA



NEO CATH[®] GS-3029

- Infant I.V cannula.
- Small flexible wings and absence of injection port permits easy handling and secure fixation.
- Provided with small projection on main body for ease of handling during insertion.
- Sterile, peel open Blister Pack.
- Box of 100/ Master Box of 1000.

NEOCARE® GS-3211

- Paediatric intravenous cannula.
- Suitable for subjects with capillary brittleness,
 where non-traumatic cannulation is a requirement.
- Provided with wing cover that assists during the insertion of cannula.
- Sterile, peel open blister pack, with tyvek Lid.
- Box of 50/ Master Box of 1000.

SPECIFICATIONS Size Length Colour 24 G 19 mm Yellow

MICRON® GS-3031

- Specially designed neonatal I.V cannula, for peripheral intravenous infusion.
- Specially designed clips for better maneuvering.
- Fine and smooth design for use in neonates to deal with the specific problems associated with small and soft vein cannulation.
- Main body is specially designed with projection for low point gap to facilitate better insertion.
- Small flexible wings and absence of injection port permits easy handling and secure fixation.
- Transparent flashback chamber for blood visualization.
- · Sterile, peel open Blister Pack.
- Box of 100/ Master Box of 1000.

SPECIFICATIONS		
Size	Length	Colour
26 G	19 mm	Violet

INFLOW SET VENTED INFUSION SET




INFLOW SET® GS-3326

- Infusion set with built in air vent and bacteria barrier hydrophobic filter.
- Sharp piercing spike for easy insertion in I.V container.
- Specially designed transparent cylindrical collapsible drip chamber to visualize the flow rate.
- Disc type 15µm fluid filter at the bottom of drip chamber filters any particulate matter in the I.V fluid.
- Provided with latex free "Y" type injection port made up of poly isoprene for additional medication.
- Efficient roller controller ensures better flow control.
- Provided with super smooth kink resistant tube for unobstructed flow.
- Fixed Luer lock for secure connection to all standard devices.
- · Sterile, individually ribbon packed.
- · Also available in blister pack.
- Polybag of 25/ Master Box of 500.

SPECIFICATIONS	
Tube Length	Drop Rate
150 cm	20 drops per ml

INTRAFLOW AS SAFETY INFUSION SET

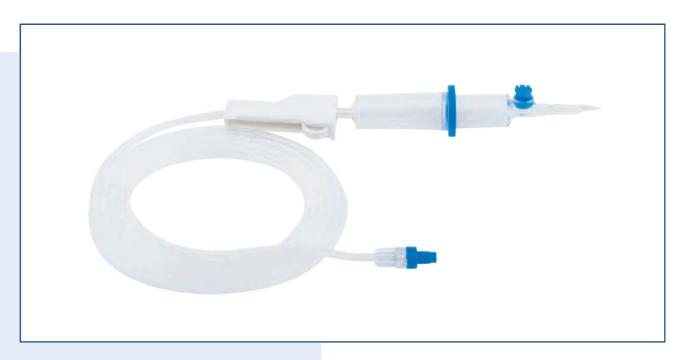


INTRAFLOW AS TM GS-3061 AS

- · Safety vented infusion set.
- Longer drip chamber fits firmly in the user's hands, ensuring a better grip.
- Unique air stop feature, which is a hydrophilic filter membrane present at the bottom of the drip chamber prevents the entry of air in case of emptied fluid container.
- Innovative prime stop, which is a hydrophobic membrane placed on the fixed luer lock fitting that stops fluid dripping on the floor.
- · Soft, clear 180 cm long kink resistant tubing.
- Modified roller controller, provided with tube docking port to hold the set when not in use.
- Spike dock is provided in roller controller for safe disposal of set after use.
- · Sterile, individually blister packed.
- · Polybag of 10/ Master Box of 200.

- ★ Air Stop Prevents Air Embolism
- ★ Prime Stop Automates Priming ★ Luer Lock

SPECIFICATIONS	
Tube Length	Drop Rate
180 cm	20 drops per ml

ROTAFLO EXTENSION SET

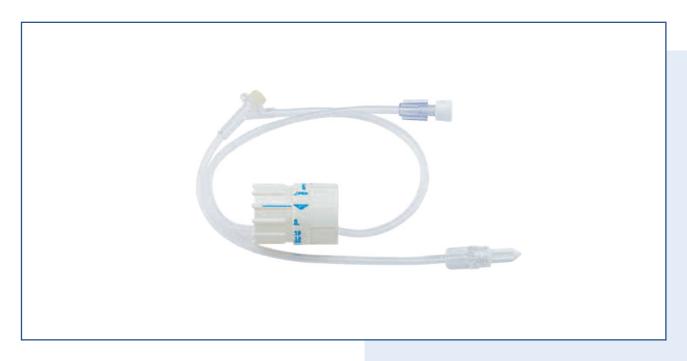

ROTAFLOTM GS-3303


- Flow regulator with extension set for precise flow control.
- Reliable flow rate accuracy across the full adjustable scale.
- Male luer connector at one end, female luer connector at other end.
- Ergonomic dial design with ridges for easy grip and adjustment.
- Precise flow rate control of IV fluids with range of 5 ml/hr-250 ml/hr.
- · Soft and kink resistant tubing.
- Provided with a latex free "Y" type injection port made up of isoprene intermittent medication.
- Can be used with all standard I.V sets using luer lock connections.
- · Sterile, individually packed.
- Box of 10/ Master Box of 100.

SPECIAL FEATURES

- ★ Flow Control Extension ★ Ergonomic Dial Design
- ★ Compatible with Standard Sets ★ Y Injection Port

Tube Length I.D/O.D 170 cm 2.9 / 4.1 mm

VOLUMETRIC PLUS/ VOLUFIX PLUS M.V SET

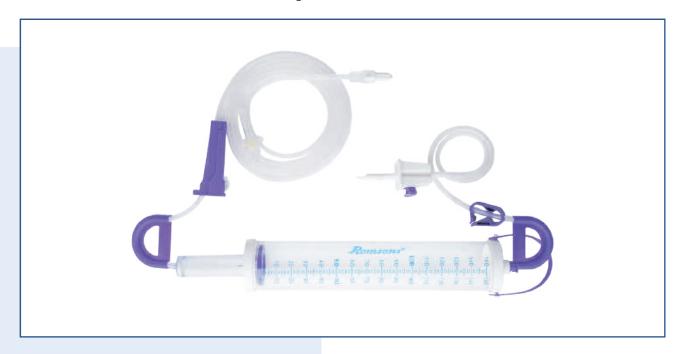

VOLUMETRIC PLUS® GS-3094 P

- · Measured volume set with burette chamber.
- Soft, large caliberated, flexible burette chamber manufactured from PETG material, suitable for most of the pharmaceutic fluids.
- Sharp piercing vented spike provided with built-in air inlet with bacterial barrier filter.
- · Micro drip with reduced drop size of 60 drops/ml.
- Easy to read scale with white contrast background.
- Floating auto shut-off valve acts as volume indicator and automatically shuts off the path to prevent air embolism.
- "Y" injection site for provision of additional medication.
- · Patented, no-kink device prevents kinking of tube.
- Sterile, individually packed.
- Box of 10/Master Box of 100.

VOLUFIX PLUS® GS-3082 P

- 100 ml burrette type Chamber with 10 ml overflow limit.
- · Micro drip with drop size of 60 drops per ml.
- Burrette chamber is made of PETG material.
 Patented design hanger facilitates hanging on I.V stand.
- Floating auto shut-off valve act as floating indicator and automatically closes the drain path to prevent air in line.
- No-kink device prevents kinking of tube during transit
- · Separate ports for extra medication and airvent.
- · "Y" injection site for additional medication.
- Roller flow controller provides precise flow control.
- Sterile, individually packed.
- Box of 10/Master Box of 100.

SPECIFICATIONS		
Capacity	Drop Rate	
150 ml	60 drops per ml	

SPECIFICATIONS	
Capacity	Drop Rate
110 ml	60 drops per ml

BI-VALVE THREE WAY STOP COCK



BI-VALVE® GS-3040

- Specially designed lipid resistant 3-way stop cock with smooth tap rotation suitable during pressure infusion and invasive blood pressure monitoring.
- · Made from clear, transparent polycarbonate.
- Rigid construction with circular flow channel prevents air entrapment.
- Minimum residual volume and uniform bore ensures maximum flow and minimum turbulence.
- Provided with two female luer lock ports and one male luer lock port.
- Screw luer lock ensures safe and secure connection to other standard luers.
- Sterile, individually packed in blister pack.
- Box of 50/ Master Box of 1000.

EXXACTA® GS-3146

- Specially designed three way stop cock with smooth rotation, suitable during pressure infusion and invasive blood pressure monitoring.
- Compatible with lipid/emulsions/ most pharmaceutics.
- · Made up of clear, transparent polymer.
- Rigid construction with circular flow channel prevents air entrapment.
- Universal 6% male/ female luer taper facilitates leak proof connections with corresponding devices.
- Male luer lock with ergonomically designed rotating nut prevents accidental disconnection of the line.
- Residual volume, not more than 0.20 ml, sustains line pressure upto 5 bar.
- Sterile ; Box of 50/ Master Box of 1000.

SPECIFICATIONS

Residual Volume	Sustainable Pressure
≤ 0.20 ml	5 Bar

- ★ Minimum Prime Volume
- ★ Lipid/ Majority Pharmaceutic Compatible

VEIN-O-LINE/ MARVEL EXTENSION LINE



VEIN-O-LINE® GS-3034

- Low pressure large bore extension line.
- Male luer lock at one end and 3-way stop cock at other end and fits all standard equipment.
- · Sterile, individually packed.

VEIN-O-LINE® L GS-3032

 Low Pressure / Large Bore Extension Line with male and female luer lock at either ends.

VEIN-O-LINE® T GS-3034 T

- Extension Line with 'T' type Luer lock connector.
- Intermittent injection port at one end, line clamp and a female luer lock connector at other end.
- · Length: 13cm; Box of 50/Master box of 500.

VEIN-O-LINE® MARVEL GS-3144

- Low pressure /large bore extension line with three way stop cock comprising one swabable needle free port for safe access.
- Swabable ports prevent thrombotic occlusion in I.V. line that inhibits microbial contamination and colonization.
- Male luer lock at one end and 3-way stop cock at other enabling the user to run two lines simultaneously.
- Sterile, individually packed in blister pack.
- Available Length: 10 cm.
- Box of 50/ Master Box of 1000.

SPECIFICATIONS

Length	Box/ Master Box
10, 50, 100, 150 cm	50/500
200 cm	20/200

- ★ 3 Way Stop Cock ★ Swabable Needle Free Port
- ★ Inhibits Microbial Contamination

PM-O-LINE/ HI LINE EXTENSION LINES



PM-O-LINE® GS-3036

- Small bore high pressure extension line.
- Suitable for high pressure monitoring and for connection between syringe infusion pump and patient.
- Multi-purpose extension line provided with male luer lock at one end and female luer lock at other end that fits all standard equipment, provides secure fitment.
- · Can sustain a pressure upto 10 kg/cm² or 140 psi.
- · Sterile, individually packed in blister pack.
- Box of 50/ Master Box of 500.

SPECIFICATIONS Length Prime Volume 50 cm 0.45 ml 100 cm 0.90 ml 150 cm 1.40 ml 200 cm 1.80 ml

HI LINE® GS-3037

- Polyethylene high pressure monitoring extension line.
- Non-reactive and highly compatible to all the commonly infused fluids.
- Small bore high pressure extension line suitable for high pressure applications.
- · Luer lock connectors provide secure fitting.
- · Can sustain a pressure upto 40 kg/cm² or 560 psi.
- Sterile, individually packed in blister pack.
- · Length: 100, 150, 200 cm
- Box of 20/ Master Box of 200.

HI LINE® SPIRAL GS-3037 S

- Specially designed spiral type extension line for greater patient comfort & freedom of movement.
- It expands and contracts gently with patient movement and protects the patient, cannula, infusion equipment from accidental damage.
- · Length: 200 cm
- Box of 20/ Master Box of 200.

EXTEENA FAMILY NEEDLE FREE DEVICES



NEEDLE FREE DEVICES

- Needleless connectors offer a safe way to access I.V lines without the use of a needle.
- Prevents thrombotic occlusion and inhibits microbial colonization.
- Secure fitment to I.V cannula, I.V set, extension lines and CVC.
- · Closed system controls infection.
- Swabable port prevents infection and allows safe access.
- · Clear PC body, Latex free.
- · Low reverse pressure facilitates easy connection.
- Compatible with MRI.
- · Sterile, blister pack.
- Box of 20/ Master Box of 200.

I.D/O.D Length 1.40/2.40 mm 15 cm

EXTEENA® UNO GS-3048U

- Needle free I.V connector line at one end and rotating luer slip lock at the other end.
- Single extension line.
- · Low prime volume: 0.30 ml.

EXTEENA® DUO GS-3048D

- Needle free I.V connector with dual extension for running two lines simultaneously.
- · Low prime volume: 0.60 ml.

EXTEENA® TRIO GS-3048T

- Needle free connector with tri extension line, rotating luer slip lock and 3 clamps for clamping.
- · Low prime volume: 0.90 ml.

EXTEENA® QUADRO GS-3048Q

- Needle free connector with four extension line.
- Low prime volume: 1.20 ml.

SAFSITE/ HEP CAP CONNECTORS



SAFSITE® GS-3018


- Needle free I.V Connector with split septum technology.
- Luer lock for secure fitment on one end and needle free port on the other.
- · Low reverse pressure facilitates easy connection.
- Swabable port prevents infection and allows safe access.
- · Clear ABS body.
- · Latex free, compatible with MRI.
- · Sterile, indivodually packed in blister pack.
- · Low prime volume: 0.08 ml.
- Box of 100/ Master Box of 2000.

HEP CAP GS-3114

- · Latex free injection stopper.
- Provided with male luer lock for snug fitment to an I.V cannula and central venous catheter.
- Integrated one way valve allows drug administration without unscrewing it.
- Swabable port prevents microbial infections and allows safe access.
- · Sterile, individually packed.
- Master Box of 2000.

- ★ Needle Free I.V Connector ★ Swabable Port
- ★ Minimizes Microbial Contaminations

CENTRO CENTRAL VENOUS CATHETER KIT

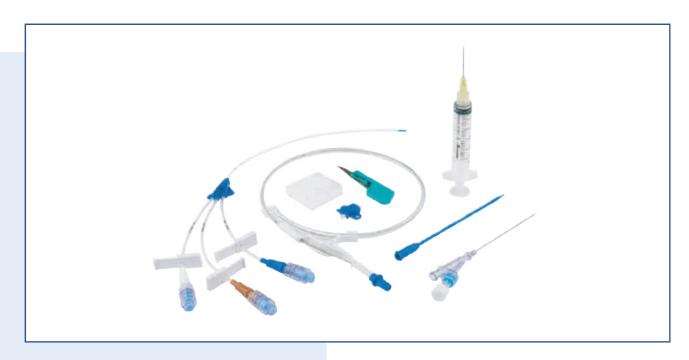

CENTRO® GS-3042

- Manufactured from specially formulated and biocompatible PU material.
- Polyurethane anti-kink soft tip ensures firm insertion and minimizes vessel erosion.
- · Radiopaque catheter with clear marking.
- Soft and flexible "J" tip guide wire prevents perforation and provides good torque to ensure firm non-traumatic insertion.
- · Based upon seldinger technique.
- Sterile, individually packed;
- · Box of 10/ Master Box of 20.

KIT CONSISTS OF:

- Soft tip Polyurethane radio-opaque catheter with luer lock connector
- · "J" tip Guide Wire with dispenser
- "Y" shaped Introducer Needle, Tissue dilator
- Suture wing, Clamp fastener, Scalpel (short)
- Extension line clamps, Swabable needle free connector
- Syringe 5ml, Hypodermic needle (20G x 25mm)
- Needle tip secure sponge

SPECIFICATIONS	
Туре	Size
Single Lumen	Fr 5 x 20 cm
Double Lumen	Fr 5 x 8 cm
Double Lumen	Fr 7 x 16, 20 cm
Triple Lumen	Fr 7 x 16, 20 cm






Diseases affecting the gastrointestinal tract, which include all the organs starting from the mouth to anus, along the alimentary canal, are the focus areas of this specialty. Our gastroenterology range of products comprises of nasogastric tube, infant feeding tube, nasojejunal tube and specialised feeding bag.





NASOGASTRIC TUBES ENTERAL FEEDING TUBES

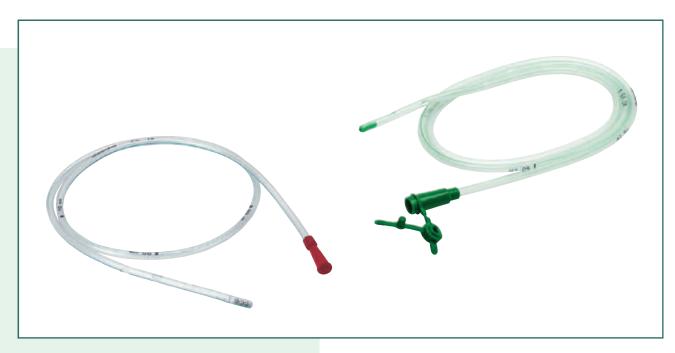

RYLE'S TUBE GS-4004

- Suitable for nasogastric introduction, for nutritional purposes or aspiration of intestinal secretions.
- Distal end is coned with corrosion resistant stainless steel balls sealed into the tube.
- · Manufactured from medical grade Polymer.
- Radio opaque line is provided throughout the length.
- Marked at 50, 60, 70 cm.
- · Available Size: 8 20 FG.
- · Sterile, individually packed.

LEVIN'S TUBE GS-4006

- · Distal end is coned with four lateral eyes.
- · Marked at 45, 55, 65, 75 cm.


SPECIFICATIONS		
Product Code	Effective Length	Total Length
GS-4004	102 cm	105 cm
GS-4032	102 cm	105 cm
GS-4006	119 cm	122 cm
	:	

ROMOLENE® GS-4032

- Nasogastric tube with coned distal end provided with radio opaque material for X-ray visualization.
- Super smooth, low friction surface facilitates smooth intubation.
- Proximal end is provided with dual type universal funnel shape connector and female luer mount with closure cap.
- Marked at 50, 60, 70 cm.
- · Available Size: 8 20 FG.
- · Sterile, individually packed.

SPECIFICATIONS		
Size	Colour Code	Box/ Master Box
6 FG	Sea Green	50/500
8 FG	Blue	50/500
10 FG	Black	50/500
12 FG	White	50/500
14 FG	Green	50/500
16 FG	: Orange	50/400
18 FG	Red	50/400
20 FG	Yellow	50/400

SALEM SUMP TUBE DUODENAL DUAL FLOW TUBE

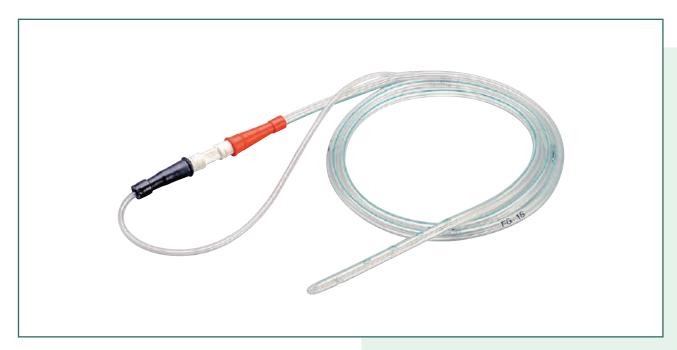


SALEM SUMP TUBE GS-4018

- Dual flow type tube is designed to prevent the accumulation of gastric secretions and air, thus allowing decompression of the stomach by continuous aspiration without trauma.
- Double lumen tube throughout the length provides air vent to break the vacuum when the stomach has been evacuated of accumulated liquid, air or small particulate matter.
- Manufactured from non-toxic, medical grade polymer.
- · Distal end is coned with ten lateral eyes.
- Tube marked at 45, 55, 65 and 75 cm from the tip.
- · Radio opaque line provided throughout the length.
- Proximal end is provided with universal funnel shaped connector along with double tapered connector for easy extension.
- The secondary lumen is provided with universal funnel shaped connector to accommodate catheter mount syringe to facilitate irrigation.
- Sterile, peelable pack.
- Box of 20/Master box of 400.

SPECIFICATIONS

Size	Effective Length	Total Length
10 FG	87 cm	90 cm
12 FG	119 cm	122 cm
14 FG	119 cm	122 cm
16 FG	119 cm	122 cm
18 FG	119 cm	122 cm
20 FG	119 cm	122 cm

INFANT FEEDING TUBE NASOGASTRIC TUBE

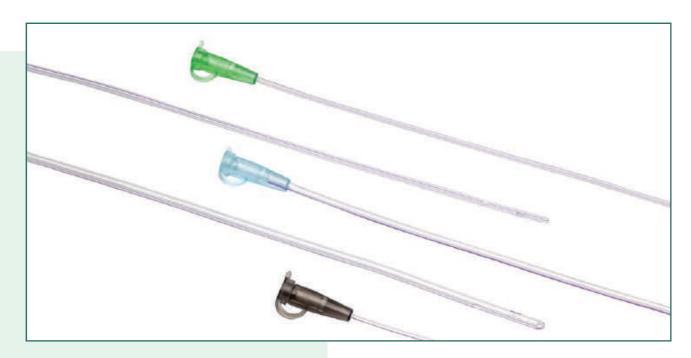

INFANT FEEDING TUBE GS-4008

- · Nasogastric, Infant feeding tube.
- The tube features low friction, super smooth frozen surface ensuring non-traumatic intubation.
- The tube is marked at 20 cm starting from tip to enable accurate placement.
- Distal end of the tube is coned for smooth intubation and equipped with two lateral eyes for efficient drainage and administration.
- Proximal end is provided with female luer mount for easy connection to feeding funnel or syringe.
- Radio-opaque line facilitates accurate tube placement.
- · Color coded for immediate size identification.
- Effective Length: 50 cm, Total Length: 52 cm.
- · Sterile, straight packed in peelable soft blister pack.
- · Box of 100/ Master Box of 2000.

SPECIFICATIONS	
Size	Colour Code
4 FG	Clear
5 FG	Grey
6 FG	Green
7 FG	Clear
8 FG	Blue
9 FG	Orange
10 FG	Black





Single-use devices are intended for use on one patient during the procedure. Not intended to be reprocessed (cleaned, disinfected/sterilized). They should be disposed of after one use, in accordance with local waste management system regulations.

RESPIROMETER RESPIRATORY EXERCISER



RESPIROMETER® GS-6018

- · Breathing exerciser.
- Innovative system to exercise respiration through inspiration.
- Specially designed with 3 stage chambers to provide for efficient exercise to patient on a step up basis.
- Device comprises of base plate with central part containing volume designated three chambers with a colored sphere.
- · Encourages deep and prolonged inspiration.
- Three spheres rise proportionately to the inspiration exerted by the user.
- Innovative design, can be dismantled into parts for cleaning and disinfection.
- Box of 1/ Master Box of 50.

UNIT CONSISTS OF:

- Main body comprising of three volume chambers
- Base plate
- · Spheres of Red, Yellow and Blue colour
- Connecting tube
- · Patented mouth piece with in-built filter

Volume
600 cc/sec
900 cc/sec
1200 cc/sec

BREAZER VOLUMETRIC RESPIRATORY EXERCISER



BREAZER® GS-6083

- Designed to encourage deep breathing exercises to restore natural breathing rhythm.
- Allows visual progress towards improving inspiratory volume.
- Provides measurement of air flow and air volume of the breathing process
- · Low resistance mouthpiece with inbuilt filter.
- · Goal indicator to indicate achieved result.
- · Made of anti-shock material with a built-in handle.
- · Latex-free, individually packed in a box.
- Box of 1/ Master Box of 20.
- UNIT CONSISTS OF:
- Transparent main body comprising of two chambers
- Connecting Tube
- · Mouthpiece with inbuilt filter


SPECIAL FEATURES

- ★ Allows View Progress ★ Goal Indicator
- ★ Encourages Deep Breathing Exercise

SPECIFICATIONS

Product Model	Volume Measurement
Breazer 2500	Upto 2500 ml

PERITONEAL DIALYSIS SET CATHETER SET

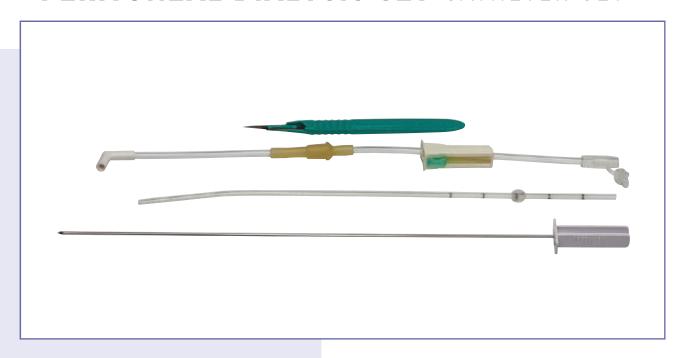

PERITONEAL DIALYSIS CATHETER SET GS-3010

- Stainless steel trocar is provided with triple angle sharpness to facilitate smooth penetration.
- Spare scalpel blade is supplied to provide the incision.
- Multihole, open end catheter is manufactured from medical grade polyamide.
- The trocar and the tip of the catheter is perfectly matched to facilitate trauma free entry.
- Latex free flash ball type injection port is provided on the junction unit for additional medication.
- Complete kit is placed in thermoformed tray, packed in peelable pouch.
- · Sterile, ready for use.
- Box of 25/Master box of 200.

UNIT CONSISTS OF:

- Stainless steel trocar with multihole smooth tip catheter.
- · Junction unit with control clamp and scalpel.

PERITONEAL DIALYSIS TRANSFUSION SET GS-3012

- Specially designed "Y" shaped transfusion set for administering dialysis solution.
- Provided with two upper control clamps to facilitate alternative change of fluid bottles.
- Two lower control clamps help easy input and drainage of solution.
- Proximal end is provided with suitable adapter for easy connection to catheter unit.
- Sterile, individually packed in peelable pouch.
- Box of 10/ Master box of 100.

SPECIFICATIONS		
Size	Perforation	Tip Marking
Adult	100 mm	190, 210, 230, 250, 270 mm
Child	50 mm	130, 150, 170mm

E.C.G ELECTRODES DISPOSABLE ELECTRODES



E.C.G ELECTRODES GS-9045

- Disposable ECG electrodes, white foam type, centred snap with conductive solid gel column ensures quality transmission.
- The electrode is designed for conditions where adhesion is difficult such as holter patient monitoring, where excessive moisture may be a problem.
- ECG solid gel electrodes eliminate the mess of liquid gels, while still offering low impedance.
- Our proprietary solid gel holds aqueous content in a polymerized hydrogel formulation, maintaining good body contact and ease of fixation.
- Non-irritating, non-sensitizing and non-toxic.
 Aggressive medical grade adhesive.
- Increased foam coating per square mm of adhesive, maximizes adhesion under diaphoretic and stress conditions.
- Hypoallergenic, latex free, non-irritating, disposable electrodes meets and exceeds all AAMI performance specifications.
- · Pack of 60/ Box of 600/ Master Box of 7200.

SPECIFICATIONS

Size	HxW
Standard	45 x 42 mm

TRIMMER PLUS PREP RAZOR



TRIMMER® PLUS GS-9024 P

- Single edge comb style prep razor intended for pre-operational procedure.
- Stainless steel blade with platinum edge and TEFLON coating.
- Steady and firm grip, contoured handle enables smooth and clean shave procedure, without cuts or irritation.
- Unique comb style facilitates non-clogging of blade due to hair or foam.
- · Meets high standards of safety and hygiene.
- Unique design helps to conveniently snap the razor from handle by applying pressure to prevent accidental re-use and proper disposal.
- Cost-effective, only one razor required for the procedure.
- · Elegant and secure pack in soft blister.
- Sterile; Box of 50/ Master box of 1200.
- · Also available in Non-Sterile version.

SPECIFICATIONS	
Dimensions	Length
Handle & Body	60 mm

CLAMPS CLAMPING DEVICES



KLIK CLAMP® GS-5084

- Suitable for clamping the umbilical cord of newly born baby immediately after their birth.
- Provided with dual security lock with distinct "click" sound to indicate the correct locking of umbilical cord.
- Finger grip ensures safe and convenient handling of the clamp, particularly when the gloves are wet.
- Provided with grooves all along the length to prevent slipping of umbilical cord and retain it in exact position.
- Manufactured from non-toxic medical grade polymer.
- · Sterile, individually packed in blister pack.
- Box of 100/ Master box of 4000.

CATH CLAMP GS-1014

- Cath Clamp is suitable for clamping any catheter or tube made from flexible polymer between
 2mm to 8mm diameter.
- Provided with distinct lock "click" which is easy to open when needed.
- Finger grip ensures safe and convenient handling, particularly when gloves are wet.
- Manufactured from non toxic medical grade polymer.
- · Supplied in green colour.
- · Other colours are available on demand.
- · Non-Sterile; Box of 100/ Master box of 2000.

- ★ Dual Security Lock ★ Distinct Click Sound
- ★ Grooves all along length ensures secure clamping

DIGNITY ADULT DIAPERS



DIGNITY® PREMIUM ADULT DIAPERS GS-8405

- · Open -Style Unisex Diapers with adjustable tapes.
- Advanced Leakage Protection for upto 10-12 Hours.
- · Suitable to be used during the day as well as night.
- · Super Absorbent Core Locks Moisture quickly.
- ADL, Acquisition Distribution Layer helps in rapid fluid dispersion.
- Leak Guard with Barrier Leg-Cuffs provides better urine and faecal containment.
- · PE Back Sheet ensures cloth-like comfort.
- Wetness Indicator indicates the need for a diaper change.
- · Pack of 10/ Master Box of 120.

DIGNITY® PREMIUM PULL-UP DIAPERS GS-8423

- Advanced Leakage Protection for upto 10-12 Hours.
- Suitable to be used during the day as well as night.
- · Easy to wear and take off.
- Stretchable Elastic Waistband.
- Super Absorbent Core Locks Moisture quickly, leaving the wearer dry and comfortable for hours.
- Inner Leg Cuffs with Leak Guard provides better urine and faecal containment.
- Breathable textile Back Sheet feels cloth-like soft on the skin.
- Wear for a contoured, discreet, and snug-fit experience.

SPECIFICATIONS	
Size	Waist
М	28 - 45 in
L	38 - 54 in
XL	48 - 57 in

SPECIFICATIONS	
Size	Waist
M/ L	24 - 55 in
L/ XL	34 - 61 in
XL/ XXL	45 - 68 in

DIGNITY DISPOSABLE UNDERPADS



DRYTEX UNDERPADSGS-8407

- Designed to protect linens and mattresses from any kind of fluid damage.
- Extra soft top sheet provides cloth like comfort.
- Super absorbent core and unique quilted pattern locks moisture quickly.
- Tear and slip-resistant, waterproof PE back sheet helps prevent any leakage.
- · Pack of 10/ Master Box of 120.

DIGNITY[®] MATTEY UNDERPADS GS-8414

- Disposable underpads designed to protect linens and mattresses from any kind of fluid damage.
- Suitable for bedridden patients or for people suffering from incontinence.
- Extra soft top sheet made from non-woven fabric provides cloth-like comfort and helps keep skin dry and healthy.
- Silicone release liners at the back help prevents underpad displacement
- Super absorbent core locks moisture quickly, leaving the user dry and comfortable for hours.
- Unique quilted pattern helps in even and rapid absorption of the fluid.
- · PE back sheet provides cloth-like comfort.

SPECIFICATIONS Size Pack/ Master Box 60 x 90 cm 10/120

SPECIFICATIONS

Size	Pack/ Master Box
60 x 90 cm	10/120
60 x 60 cm	20/240
60 x 40 cm	25/ 400

DIGNITY WET WIPES



SPONGEETM BODY WIPES GS-8425

- Gently Cleanses and Moisturises skin with a mild refreshing fragrance.
- Crafted from Hypoallergenic fabric, Spongee helps prevent any skin rashes or infections.
- Enriched with Aloe Vera and Vitamin E, these wipes keep skin happy and healthy.
- Clinically Proven, Alcohol & Paraben-Free, Spongee helps maintain the pH balance of the skin.
- Ideal for bedridden patients or for post operative use in hospitals, nursing homes and home care. Use it before changing diapers.
- · Microwavable ; can be used hot or cold.
- · Master Box of 48 Packs.

SPECIFICATIONS Size Pack 240 x 300 mm 10 pcs 300 x 300 mm 10 pcs

SPARKLEAN SURFACE CLEANSING WIPES GS-8455

- · Multipurpose Surface Disinfectant wipes.
- Can be used to sanitize multiple surfaces such as Gadgets, Door Knobs, Steering Wheels and so on.
- Kills 99% Bacteria & Viruses ; Nonresidual.
- · Paraben Free ; pH Balanced.
- · Master Box of 24 Packs.

KLEANZO HAND CLEANSING WIPES GS-8457

- Sanitises & Nourishes Skin.
- · Kills 99% Bacteria & Viruses.
- · Crafted from Organic Bamboo Fabric.
- Soft & Gentle ; Enriched With Aloe Vera.
- Clinically Tested; Hypoallergenic Fabric.

SPECIFICATIONS Product Code Size Pack GS-8455 240 x 300 mm 50 pcs GS-8457 150 x 200 mm 72 pcs

DIGNITY BABY WIPES



SPONGEE TH BABY WIPES GS-8428

- · Gently cleanses and moisturises the baby's skin.
- · Clinically proven; delightful mild frangrance.
- Safe to be used on baby's face and intimate areas
 Enriched with Aloe Vera & Vitamin E
- · 100% Alcohol & Paraben-Free.
- · Hypoallergenic fabric prevents skin rashes.
- Master Box of 24 Packs.

99% WATER WIPES GS-8450

- · Formulated from 99% ultra-pure water.
- · Clinically proven ; maintains the pH balance.
- Enriched with goodness of aloe vera & Vitamin E.
- · Hypoallergenic fabric ; fragrance-free.
- · 100% Alcohol & Paraben-Free.
- · Disposable, Convenient and Travel-Friendly.

SPECIFICATIONS

Product Code	Size	Pack
000.20	150 x 200 mm	72 pcs
GS-8450	150 x 200 mm	80 pcs

POOCHIE® BABY WIPES GS-8445

- Crafted from exceptionally soft 100% Biodegradable Bamboo Fabric.
- Safe and gentle on your baby's sensitive skin and the environment
- · Enriched with Aloe Vera and Vitamin E.
- Safe to be used on baby's face and intimate areas;
 these are ideal for a nappy change, hand wipe,
 and face wipe.
- Hypoallergenic fabric helps prevent any skin rashes or infections
- · Maintains the pH balance of the skin.
- · Lid pack retains moisture for a long time.
- · Disposable, Convenient & Travel-Friendly.
- Mater Box of 24 Packs.

SPECIFICATIONS

Size	Pack
150 x 200 mm	80 pcs

DISPO GUARD/ COMFIT DISPOSABLE FACE MASK



DISPO GUARD® GS-6107

- Disposable Protective Face Masks for one-time use.
- · 3-Ply Design with Melt Blown Filter.
- Fluid Repellent, filters 99% Bacteria released in droplets and aerosols.
- Extremely Soft Non-Woven Ear Loops for a comfortable wear for longer hours.
- · Adjustable Nose Bridge fits the shape of the nose.
- · Skin-friendly & breathable.
- Made from latex-free, hypoallergenic non-woven fabric.
- · Ideal for daily use.
- Box of 50/ Master Box of 2000.

COMFIT® 3D FACE MASK GS-6146

- Disposable protective face mask, ideal for one-time use only.
- 3-Dimensional construction fits all face-types perfectly, giving enough breathing space.
- 4 layer protection with 2 layers of melt-blown filter.
- Provides effective protection against bacteria released in droplets & aerosols.
- Soft & breathable, made from latex-free, hypoallergenic non-woven fabric.
- · Incredibly soft non-woven ear loops.
- · Snug-fit with complete protection.
- · Can be comfortably worn for longer time.
- Box of 25/ Master Box of 2000.

SPECIAL FEATURES

- ★ Hypoallergenic 3 Ply Mask with Melt Blown Filter
- ★ Extremely Soft Non-Woven Ear Loops

- ★ 3 Dimensional Fit ★ 4 Layer Protection
- ★ 2 Layers of Melt Blown Filter ★ Soft Ear Loops

KARE MASK/ PLUS N95 RESPIRATOR MASK



KARE MASK TM GS-6103

- Advanced double layers of melt-blown filters protects from 99% of bacteria and viruses shed in droplets and aerosols.
- · Built-in nose clip allows snug fit.
- Stretchable ear loops for comfort wear.
- · Breathable. non-woven fabric.
- · Keeps you safe from dust or pollutants.
- · Individually packed; Box of 25/ Master Box of 300.

KARE MASK TM CHILD GS-6103 C

- · Specially designed for kids.
- Light weight N95 Respirator Mask.
- Comfortable, adjustable ear loops.
- · Attractive cartoon character pattern.
- · Individually packed.

SPECIAL FEATURES


- ★ 99% Bacteria Filtration Efficiency
- ★ 6 Layer Protection ★ Soft Breathable Fabric

KARE MASK TM PLUS GS-6108

- N-95 respirator mask with double layer melt
 -blown filter.
- Adjustable Head Loops guarantee an extremely comfortable snug fit experience for longer hours.
- Ideal for surgeons and paramedical personnel and the general population.
- · Built-in Nose Clip fits the Shape of the Nose.
- Skin Friendly & Breathable. Made from Latex-free,
 Hypoallergenic Non-Woven Fabric.
- Protection from 99% bacteria, viruses, dust, and pollutants.
- Individually packed; Box of 25/ Master of 300.

- ★ 99% Bacteria Filtration ★ 6 Layer Protection
- ★ Adjustable Head Loops ★ Soft Breathable Fabric

C.V.P MANOMETER C.V.P MEASURING DEVICE

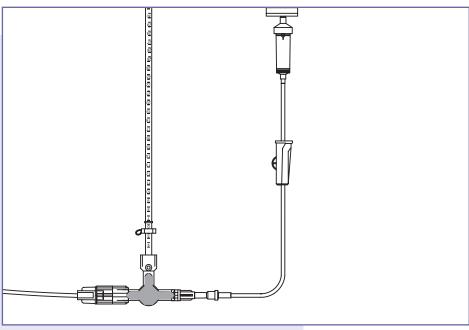

C.V.P MANOMETER® GS-3008

- Suitable for continuous or intermittent monitoring of central venous pressure during infusion, particularly in ICU, CCU, OT's and casualty departments.
- Manometer tube with graduation from -4 cm to
 +34 cm and attached to three way tap.
- Provided with 120 cm long extension tube, fitted with female luer lock on one end and male luer lock on the other end for safe connection to the central venous catheter.
- Latex free injection site is provided for extra medication.
- Sliding indicator provided to record previous reading.
- Moulded clamps with latex free rubber strings for easy fixation on the intravenous stand.
- · Sterile, individually packed in peelable pouch.
- · Individually packed in a box.
- Box of 1/ Master Box of 50.

