

SEESHEEN

see the endo-sheen

Edition 01/2018



Flexible Video Endoscopes
PRODUCT CATALOGUE



see the endo-sheen
开启可视医疗新境界

珠海视新医用科技有限公司

Mutual respect, trust and benefit

Vision

- Make the medical endoscopy solutions to be more accessible and affordable

Mission

- To be a leading manufacturer in slim flexible CMOS videoscope field





CMOS CHIP-ON-TIP

No fibers



Durable

No fiber design
More flexibility
Less maintenance



Slim

The slimmest outer
diameter reaches
2.8mm



Portable

Portable design
for anytime and
anywhere



Compact

Light source
Display & Processor
all-in-one

attach **water proof cap**
to achieve **IPX7**

for fully immersing disinfection



Completed


Tracheal Flexible Intubation Solutions

- Slim, standard and therapy type for optimized usability;
- Chip-on-Tip video technology;
- Portable design is applicable to intubation in anywhere;
- Photograph and video recording;
- Standard configure 3inch attachable display;



SPECIFICATIONS

Model	LF-X10	LF-XB10	LF-20	LF-XB20
View Angle	90°	90°	90°	90°
Insertion Tube OD	3.1mm	3.9mm	5.2mm	5.8mm
Channle ID	1.2mm	1.5mm	2.2mm	2.6mm
Working Length	610mm	610mm	610mm	610mm



Anesthesiology



Occurrences

Safe and fast Intubation

Seesheen Intubation Scope standard attaches removable 3 inch display which offers open vision and brilliant user experience. The product line offers ideal and completed intubation solutions by the features of versatile, efficient and scatheless. They can be used for the situations:

- Difficult airway intubation
- Special patient postures
- Awaking intubation
- DLT/BB position checking
- Nasotracheal Intubation
- Teaching and practice

Bedside

Flexible Bronchoscopy Solutions

- Slimmest 3.2mm outer diameter;
- Chip-on-Tip video technology;
- Left & right 120 °insertion rotation;
- Photograph and video recording;
- Portable design offers more convenience for routine examination and therapy;



SPECIFICATIONS

Model	VB-332	VB-349	VB-359
View Angle	120°	120°	120°
Insertion Tube OD	3.2mm	4.9mm	5.9mm
Channle ID	1.2mm	2.2mm	2.8mm
Working Length	610mm	610mm	610mm



Pulmonary



Occurrences

Make bronchoscopy become more accessible

Nowadays, Bronchoscopy is becoming a crucial tool of lung disease Interventional therapy. Seesheen flexible video bronchoscopes including 3 types which offer more accessible bronchoscopy therapy solutions.

- Interventional laser/cryo-therapy;
- Remove of Mucus / foreign bodies;
- Placing airway stent / Bronchial dilatation
- Lavage such as BAL/WLL;
- Dose purpose for instance fistula blocking or EBTB

E.N.T.

Flexible Rhinolaryngoscopy Solutions

- Portable design for OPD daily using;
- Wide 120° viewing angle;
- Chip-on-Tip video technology;
- Integrated LED light source;
- 130° up and down bending angles;
- Less stimulus and more accurate for examination and therapy;
- Photograph and video recording;



SPECIFICATIONS

Model	EV-439	EV-349
View Angle	120°	120°
Insertion Tube OD	3.9mm	4.9mm
Channel ID	-	2.2mm
Working Length	320mm	365mm



E.N.T



Occurrences

Softer and safer

Flexible rhinolaryngoscopy is currently one of the most commonly used techniques for visualizing the nasal cavity, the sinuses, the pharynx and the larynx. Seesheen flexible video Rhinolaryngoscope EV-349 which has 2.2mm working channel can be widely used for:

- Polypectomy of nasal and vocal core
- Adenoids removal
- Rhinopharyngocele checking and removal
- Rhinopharyngocele check and removal
- Accurate dosing
- Checking bleeding points

RIRS & PCNL

Flexible Endoscopy Solutions

- Better vision and durability than traditional fiber scopes;
- Video Chip-on-Tip technology;
- 275° bending to reach lower calyx;
- Left & right 120° insertion rotation;
- Photograph and video recording;
- Video signal output ability to medical grade display;



UR-332




CN-355



SPECIFICATIONS

Name	Cysto-nephroscope	Uretero-renoscope
Model	CN-355	UR-332
Insertion Tube OD	16.5F	9.6F
Channel ID	7.8F	3.6F
Working Length	380mm	670mm

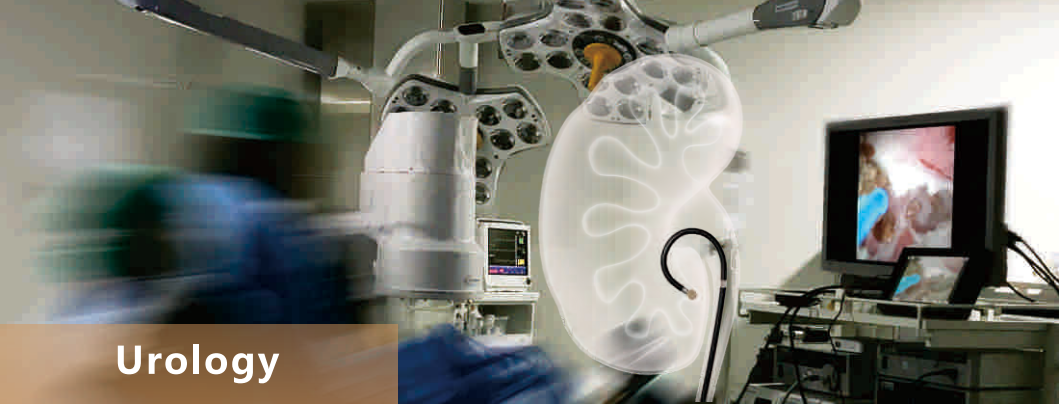

Urology



Occurrences

RIRS

Retrograde Intra renal surgery

- RIRS is an effective surgery which is using flexible uretero-roscope to enter kidney retrogradely. Seesheen UR-332 flexible uretero-roscope combines chip-on-tip technology which gives better vision and durability than the fiber scopes.

PCNL

Percutaneous Nephrolithotomy

- Compared with rigid cysto-scope, Seesheen CN-355 offers soft and flexible cystoscopy for examination and therapy; it also can be used for PCNL when the rigid scopes meet blind area during the operation.

Flexible video choledochoscope

LTCBDEC

Choledochoscopy Solutions

- Video Chip-on-Tip technology;
- Photograph and video recording;
- Left & right 120° insertion rotation;
- Video signal output ability to medical grade display;



SPECIFICATIONS

Model	CV-332	CV-349
View Angle	120°	120°
Insertion Tube OD	3.2mm	4.9mm
Channel ID	1.2mm	2.2mm
Working Length	450mm	380mm



Occurrences

LTCBDEC

Laparoscopic transcystic common bile duct exploration with choledochoscope

- Coordinate with laparoscope, Seesheen choledochoscope can be used for LTCBDEC which is an ideal solution for lithotripsy. CV-332 with 3.2mm outer diameter and 1.2mm instrument channel can support doctor to search and break stones in smaller biliary duct.



CE 0197



SEESHEEN

see the endo-sheen



ZHUHAI SEESHEEN MEDICAL TECHNOLOGY CO., LTD.

ADD: 4f, bldg. b, 6th plant, Pingxi 10th road, Nanping science & technology industrial park, Zhuhai, 519060, China.
TEL: +86-0756-8676359 FAX: +86-0756-8676349 EMAIL: intl@seesheen.com
WEB: WWW.SEESSHEEN.COM

all rights reserved. Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.