

# ENDO ECG-301

3-Channel ECG



# ENDO ECG-301

## 3-Channel ECG


ENDO ECG-301 is the new generation of 3-channel ECG, which incorporates the benefits of the advanced technology and innovative design in an extremely compact and lightweight device. It is specifically developed for those who pursue superb portability without compromising in the quality for mobile environments.



Ultra-portable Design



Full-touch Operation

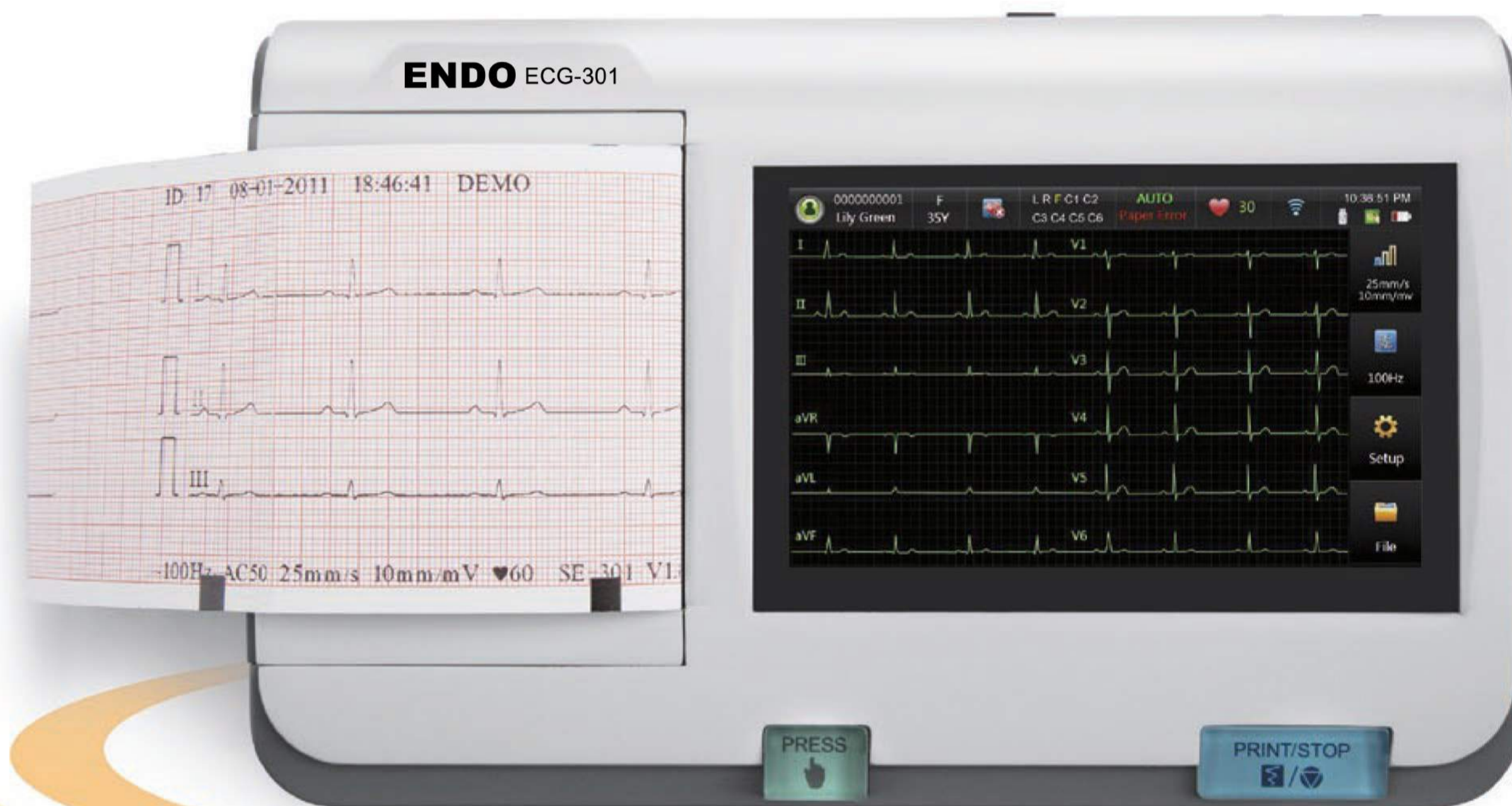

## All for a reliable result

### CMRR

A high Common Mode Rejection Ratio of more than 140dB restrains signal noises to the greatest extent, thus offer you traces with excellent quality.

### Band width

A wide band width from 0.01 to 300 Hz guarantees the accuracy of the ECG signal as well as the interpretation especially for pediatric patients.

### Sampling Rate

The leading sampling rate as high as 16,000 Hz aims at capturing all the details to generate precise result.

## Convenient Information Communication

### Powerful External Device Support

ENDO ECG-301 arms with two USB ports as well as one microSD card slot, which far extends its capability by supporting various external devices such as barcode scanner, external printer, USB Memory Drive as well as microSD cards.

### Bi-directional Communication with Data Management System

Via the LAN or build-in Wi-Fi (optional), ENDO ECG-301 is able to receive orders from ENDO ECG-301 Data Management System and transmit ECG reports to it to have further analysis, editing, saving or printing.

### Multi-format Report (Optional)

ENDO ECG-301 is able to export ECG reports in standard PDF format and optional SCP, XML and DICOM format, which offers you greater capability in integrating the ECGs into different data management systems.



### 120s waveform freeze

The freeze function enables you to review 120s waveform and choose any typical segment to analyze and print, which helps to generate integrated reports.

### Signal Quality Detection

By showing different colors to indicate normal, poor or missing signals, this function enables you to identify the signal quality even at a glance.

### Reliable SEMIP® ECG algorithm

SEMIP® ECG algorithm offers you an accurate and reliable reference to back you up when making diagnosis.



## ENDO ECG-301

Multi-channel ECG

### About ENDO

ENDO INDONESIA was founded in 2006 as a company engaged in importation and distribution of Medical Equipment. Today our products include Electromedical Equipment, Reagents, Reusable and Disposable Medical Devices, among many others and accessories. And now we had grown into the manufacturing field of Medical Recording Chart Paper.

ENDO is a company that fulfills legal requirements and receives approval from company's legal entity of Indonesia, with number C-05336 HT.01.01.TH.2006. Supported by experienced technicians and proper equipments, we also provide repairs and calibration for electromedical equipment. ENDO exist to serve your various needs of Medical Equipments and highly concerned about customer satisfaction to promote long term working relationship.



#### PT. ENDO INDONESIA

Jl. Raya Menganti No. 14, Wiyung - Surabaya 60223

Tel : (031) 767 3636 | Email : [info@endo.id](mailto:info@endo.id)

Fax : (031) 767 3737 | URL : <https://endo.id>

Toll Free : 0800 177 ENDO (3636)



Certified Company  
No. 16 00 A 14006